WARNUNG:

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**WARNING:**

We strongly advise you to read the comprehensive instructions for use including all warning and safety notes before using the device for the first time! These brief instructions do not replace the comprehensive instructions for use!

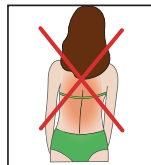
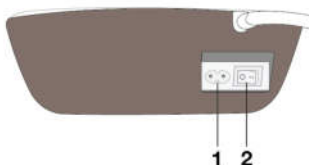
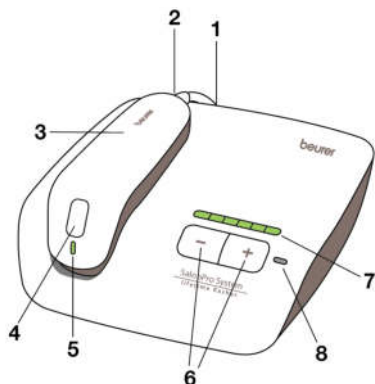
Remember that it may take a complete hair growth cycle (18-24 months) for hair removal to be permanent.

**WARNING:**

Do not treat the same area of skin more than once per hair removal session! Treating an area of skin more than once during a single hair removal session can increase the probability of side effects. Try to avoid overlapping pulses.


Do not use the device on tanned skin or if you recently suffered from sun burn. Wait two weeks following sunbathing before you use the device.

**Unit description**


- 1 Base unit socket
- 2 Power switch
- 3 Applicator
- 4 Trigger button
- 5 Indicator lamp
- 6 Energy selection buttons
- 7 Energy indicator lamp
- 8 Indicator lamp for "Cartridge empty"

Use**1**

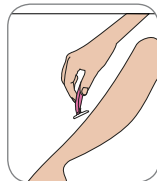
Before using the device you must shave, clean and dry off affected areas as well as keep said areas dry (do not apply perfume, make up, powder, creams, etc.).






**2**

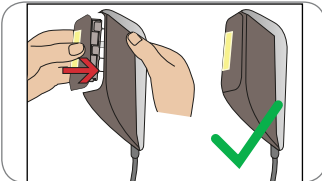Check that the cartridge is correctly plugged into the applicator.

**3**

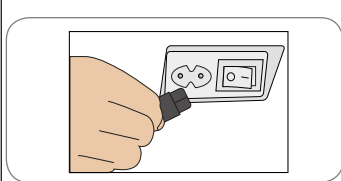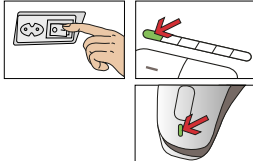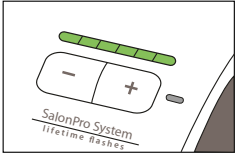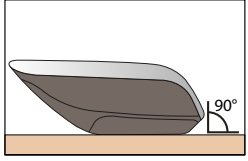
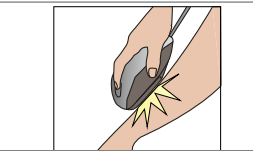
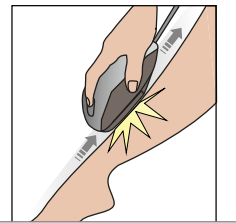
Insert the mains cable into the base unit socket. Insert the opposite end of the mains cable into the socket.

4 Switch the mains switch on the rear of the base unit -> indicator lamp and energy indicator lamp light up. 	5 Select the desired energy level using the energy selection buttons “-” and “+”. 	6 Position the light surface of the cartridge on the desired treatment surface and make sure it is fully in contact with the skin. 
7 Press the trigger button on the applicator. A bright flash appears and you will hear a brief sound. Lift the light surface of the cartridge off the skin and position it on the next treatment surface. Press the trigger button on the application as soon as the indicator lamp on the applicator lights up green. 	8 <u>Flash-glide function</u> Press and hold the trigger button = automatic flash every 3 seconds. 	

i In the case of body hair treatment with the beurer SalonPro System, the first four sessions should be performed two weeks apart. The fifth to seventh sessions inclusive can be performed four weeks apart. Any further sessions should be performed if the hair has grown back until the desired results have been achieved. In the case of treatment of hair on the face below the cheekbones, the first six sessions should be performed two weeks apart. The seventh to twelfth sessions can be performed four weeks apart. Any further sessions should be performed if the hair has grown back until the desired results have been achieved.

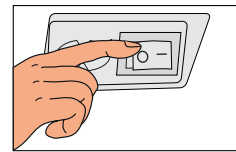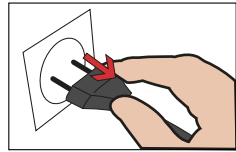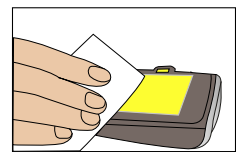After application _____

1 Switch the mains switch on the rear of the base unit -> device switched off. 	2 Remove the mains cable from the socket. 	3 Clean the light surface of the cartridge using a dry, clean cloth. 
---	--	--

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