
DE Paraffinbad Gebrauchsanweisung	4	PL Wanienska do parafiny Instrukcja obsługi	70
EN Paraffin bath Instructions for use	14	NL Paraffinebad Gebruiksaanwijzing	79
FR Bain de paraffine Mode d'emploi	23	DA Paraffinbad Betjeningsvejledning	88
ES Baño de parafina Instrucciones de uso	33	SV Paraffinbad Bruksanvisning	97
IT Bagno di paraffina Istruzioni per l'uso	42	NO Paraffinbad Bruksanvisning	106
TR Parafin banyosu Kullanım kılavuzu	51	FI Paraffiinallas Käyttöohje	115
RU Парафиновая ванночка Инструкция по применению	60		

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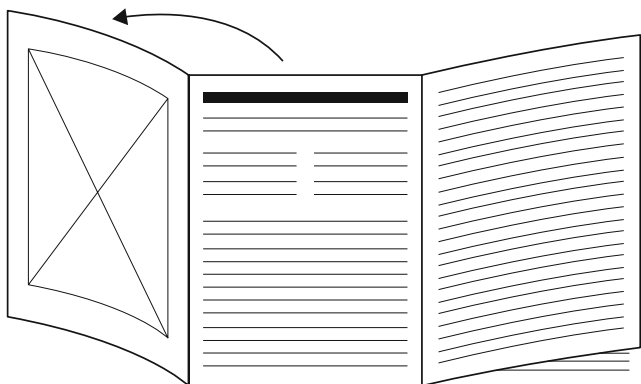Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

TABLE OF CONTENTS

1. Warnings and safety notes.....	14
2. Signs and symbols.....	16
3. Included in delivery	17
4. Device description.....	17
5. Initial use	18
6. Usage	19
7. Cleaning and maintenance.....	20
8. Accessories and replacement parts.....	21
9. What if there are problems?	21
10. Disposal.....	22
11. Technical specifications	22
12. Warranty	22

1. WARNINGS AND SAFETY NOTES

▲ WARNING

- The device is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- Do not use the device if it shows signs of damage or does not function properly. In such cases, contact Customer Services.
- Repairs must only be carried out by Customer Services or authorised retailers.

- Under no circumstances should you open or repair the device yourself, as correct functioning of the device can no longer be guaranteed should you do so. Failure to observe this will invalidate the warranty.
- The device is intended solely for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use. This device is not suitable for commercial use or for use in the medical domain. Only use the device in the private domain.
- Do not use the device if you have heart disease, flushed skin, inflammations, swelling, sensitive areas, diabetes, open wounds, burns, eczema, psoriasis or reduced sensitivity to pain.
- Do not use the device if you are not sensitive to heat.
- If you are pregnant or suffer from health problems, please consult your doctor before using the device.
- If you suffer from acute inflammations or a peripheral vascular disease, please consult your doctor before using the device.
- Never operate the paraffin bath without paraffin wax.
- Remove jewellery from your hands and feet before use.
- Never heat the paraffin wax in its plastic packaging. Never heat the paraffin wax on the hob, over a naked flame or in a microwave. Only heat the paraffin wax in the paraffin bath.
- Children do not recognise the danger posed by electrical devices. Make sure that the device cannot be used by children without supervision.
- Check that the device and the mains cable are in perfect working order before each usage. In the event of damage, the device must not be used.
- Before connecting the device to the mains, check whether the local mains voltage corresponds to the voltage specified on the device.
- Do not use the device near water (e.g. bath, wash basin etc.).
- Do not use the device while bathing, showering or sleeping.
- Never place the device in water or other liquids. Store the device in a safe place.
- Never place the device on easily damaged surfaces.
- When using the device, ensure that it is placed on a stable and level surface.
- Keep the device and mains cable away from hot surfaces. Do not place it on top of other devices, as the hot paraffin wax may cause damage.
- Do not use the device in the vicinity of aerosol products or equipment that emits oxygen.
- Do not use the device outside. The device is only to be used inside.
- Before use, ensure that the mains cable does not pose a trip hazard.
- Do not use the device in bed, whilst sleeping or if you are tired.
- Never use the device under a blanket or towel etc.
- Only use the device with the supplied accessories.
- Do not overfill the paraffin bath. Ensure that the liquid paraffin wax does not exceed the MAX fill indicator line.
- If the paraffin wax has completely melted, the device may be very hot.
- Before submerging your hands, feet or elbows in the liquid paraffin wax, use the inside of your wrist to check that the temperature of the wax is comfortable.
- Should a fault occur during use, switch off the device immediately and remove the mains plug from the socket.
- Keep the hot paraffin wax away from your eyes and other body parts that are sensitive to heat.

- If the paraffin wax permanently feels too hot or too cold against your skin, this may be a physical problem. If this is the case, consult your doctor.
- Make sure that the device is always filled with paraffin wax up to the MIN fill indicator line during use.
- Do not apply the paraffin wax to the treatment area for longer than 20 minutes.
- Do not shake or move the device during use. You could sustain scald or burn injuries.
- If there is hot paraffin wax in the paraffin bath, do not pour in any water.
- If pain, irritation, dermatitis or discomfort occurs, terminate the treatment and consult your doctor.
- Should the device accidentally fall into water, disconnect the mains plug immediately. Have the device checked by a qualified electrician before it is used again.

CAUTION

- Switch off the device after use by removing the mains plug from the socket.
- Never remove the mains plug from the socket by pulling on the cable or with damp hands.
- Do not use any sharp or pointed objects to remove the wax block from the paraffin bath. You could damage or scratch the paraffin bath.
- If you are not going to use the device for an extended period of time, store it in a cool, dry and dust-proof place.
- If you should have any questions about using our devices, please contact either your retailer or Customer Services.

INTENDED USE

This device is intended exclusively for use on the hands, feet or elbows of humans. Do not use the device on animals!

2. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	Read the instructions
	Indicates a potentially imminent danger. If it is not avoided, there is a risk of death or serious injury.
	Indicates a potentially imminent danger. If it is not avoided, slight or minor injuries may result.
	Indicates a potentially harmful situation. If it is not avoided, the system or something in its vicinity may be damaged.
	CE labelling This product satisfies the requirements of the applicable European and national directives.
	Manufacturer
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE

	Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU
	United Kingdom Conformity Assessed Mark
	Protection class II device The device is double-insulated and is, therefore, in protection class 2
	For indoor use only
	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	Importer symbol

3. INCLUDED IN DELIVERY

Check that the device packaging has not been tampered with and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

1 x MP 70 paraffin bath

1 x lid

1 x protective grille

2 x 450 g bags of paraffin wax

30 x plastic sleeves

4. DEVICE DESCRIPTION

The corresponding drawings are shown on page 3.

1 Paraffin bath

2 Lid

3 Green operation light

4 Temperature controller







5 Red heating light

6 Protective grille

7 Paraffin wax

8 Plastic sleeves

5. INITIAL USE

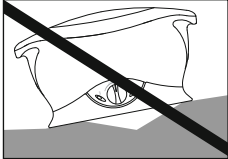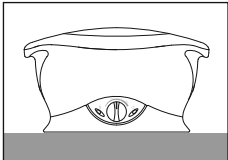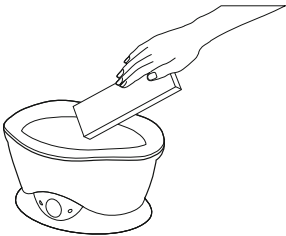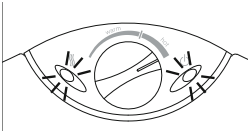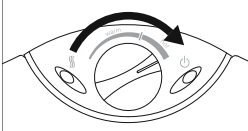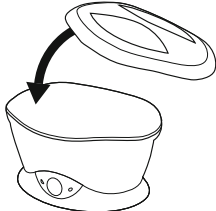• Remove the packaging material. • Make sure that all parts specified in the Included in delivery list are contained (see 1. Included in delivery). • Check the paraffin bath and mains cable for damage. • Place the paraffin bath on a flat, secure, even and water-resistant surface or table. • Make sure that the paraffin bath is empty before filling it with paraffin wax.	 
• Remove the paraffin wax from its plastic packaging. • Put the desired amount of paraffin wax (minimum 0.9 kg, maximum 1.35 kg) into the paraffin bath. In order to make it easier to put the paraffin wax into the paraffin bath, break up the blocks of paraffin wax. ⚠ CAUTION Make sure that the paraffin bath is always filled with paraffin wax up to at least the MIN fill indicator line before use.	
• Insert the mains plug into the socket (240 V~, 50/60 Hz). The green operation light on the paraffin bath will light up. The device is now switched on. The red heating light on the paraffin bath will also light up, and will only go out when the specified heat setting has been reached.	
• Turn the temperature controller to “hot” to melt the paraffin wax.	
• Now place the lid on the paraffin bath and wait until the paraffin wax has completely melted (around 2-3 hours). ⚠ WARNING If the paraffin bath is set to “hot”, never dip your hands, feet or elbows in the paraffin wax.	


- As soon as the paraffin wax has completely melted, remove the lid.
- Turn the temperature controller to the marked line between “warm” and “hot”.
- Lay the protective grille included in the delivery contents in the liquid paraffin wax and wait until it has sunk to the bottom of the paraffin bath.
- Put the lid back onto the paraffin bath and wait approximately one hour, until the paraffin wax has cooled to an ideal operating temperature (60°C).

⚠ WARNING

The room temperature can affect the power of the device and consequently the temperature of the wax.

Before submerging your hands, feet or elbows in the paraffin bath, use the inside of your wrist to check that the temperature of the wax is comfortable.

If a skin should form on the wax surface, increase the temperature by turning the temperature controller clockwise. If the paraffin wax is too hot for you, reduce the temperature by turning the temperature controller anti-clockwise.

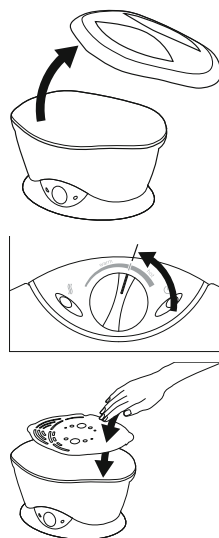6. USAGE

⚠ WARNING

Before submerging your hands, feet or elbows in the paraffin bath, use the inside of your wrist to check that the temperature of the wax is comfortable.

⚠ WARNING

Take off your jewellery (e.g. watches, rings, etc.) before treatment.


After you have completed the necessary preparations from “Chapter 6 Initial use”, you can begin using the device. To start the treatment, carry out the following steps:



- We recommend applying a lanolin-based hand lotion to the areas to be treated before treatment. This will help you achieve a better result, and will make it easier to remove the paraffin wax after treatment.

⚠ WARNING

When using the paraffin bath on your feet, only use it sitting down. Never stand in the paraffin bath.

- Briefly submerge the area to be treated (your hand, foot or elbow) for 1-2 seconds in the liquid paraffin wax. Remove your hand again and wait briefly until the wax has solidified. During this time, do not move the treated area, so that the wax can solidify properly.

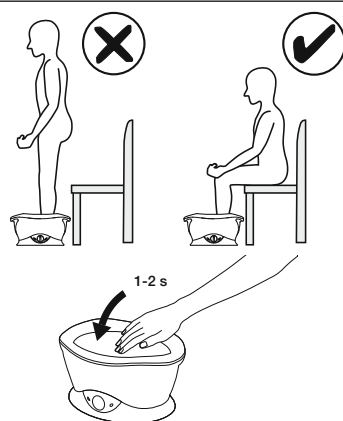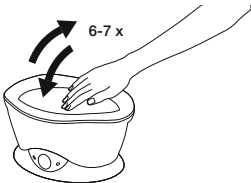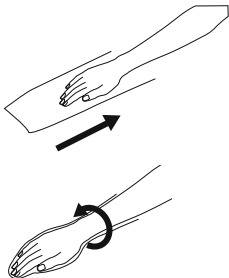
Repeat the process of submerging and pulling out six to seven times, until several layers of paraffin wax have formed over the treatment areas. If the paraffin wax is either too hot or too cold for you during the treatment period, turn the temperature controller clockwise or anti-clockwise, to increase or decrease the temperature as required.	
Now put a sheet on the treated area, press it onto the paraffin wax and twist the end shut. Keep the sheet on for 20 minutes. ⚠ WARNING If you are wearing the plastic sleeves on your feet, do not run around in them. You could slip.	
After 20 minutes, take the sheet and the layers of wax off the treated area. Dispose the sheet in the general waste; it cannot be reused.	
After each use, take the protective grille out of the paraffin bath.	

7. CLEANING AND MAINTENANCE

⚠ CAUTION

If the device is only being used by one person, change the paraffin wax at regular intervals and clean the paraffin bath whilst doing so.

If the device is being used by several people, the paraffin wax should be changed before each treatment for reasons of hygiene.

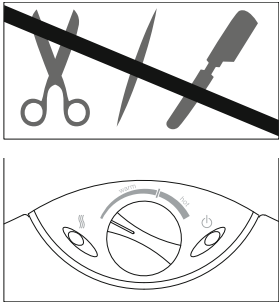- Pull the mains plug out of the socket. - Take the protective grille out of the paraffin bath. - Allow the paraffin wax to solidify overnight.	
⚠ WARNING Do not use any sharp or pointed objects to remove the wax block from the paraffin bath. You could damage or scratch the paraffin bath. If the paraffin wax has solidified, press on one end of the wax block and lift it out of the paraffin bath. If the wax block cannot be removed easily, turn the temperature controller to “warm” for ten minutes. Then remove the mains plug and take the wax block out of the paraffin bath.	

- Then wipe the inside of the paraffin bath dry with a clean towel.

⚠ CAUTION

Never fill the paraffin bath with water in order to clean it.

- After cleaning, store the paraffin bath in a cool, dry and dust-proof place.

8. ACCESSORIES AND REPLACEMENT PARTS

To purchase accessories and replacement parts, please visit www.beurer.com or contact the corresponding service address (according to the service address list) in your country. Accessories and replacement parts are also available from retailers.

Designation	Item and/or order number
2 x 450 g of paraffin wax + 30 plastic sleeves	589.43

9. WHAT IF THERE ARE PROBLEMS?

Problem	Possible cause	Solution
The device is not working		For repairs, please contact Customer Services or an authorised retailer
Heating lamp/ operation lamp does not light up	No power supply	Plug the mains adapter into the plug socket
		The mains adapter is faulty. Contact Customer Services or your retailer.
	LEDs have reached the end of their service life. LEDs faulty	For repairs, please contact Customer Services or an authorised retailer
Hard wax in paraffin bath	The oil in the paraffin wax is used up	Renew the paraffin wax. Only use genuine Beurer paraffin wax
	May occur due to frequent use	
	A paraffin wax with a higher melting point is used	
Discoloured paraffin wax	May occur due to frequent use	Renew the paraffin wax. Only use genuine Beurer paraffin wax
The paraffin wax takes longer to melt than usual	The paraffin wax is too cold	Store paraffin wax at room temperature
	The device is not working	For repairs, please contact Customer Services or an authorised retailer.

10. DISPOSAL

Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).

If you have any queries, please contact the appropriate local authorities.

11. TECHNICAL SPECIFICATIONS

Model no.	MP 70
Voltage supply	220-240 V~, 50/60 Hz
Output	80 W
Protection class	II / 
Dimensions (L x W x H)	33 x 23 x 18.5 cm
Weight	1 kg

Subject to technical changes

12. WARRANTY

Further information on the warranty and warranty conditions can be found in the warranty leaflet supplied.

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UK: Beurer UK Ltd., Suite 16, Stonecross Place, Stonecross Lane North, WA3 2SH Lowton, United Kingdom
Beurer GmbH, Söflinger Str. 218, 89077 Ulm, Germany, www.beurer.com