

<b>DE</b> Shiatsu Fußbad Gebrauchsanweisung .....	4	<b>PL</b> Wanienka do wodnego masażu stóp Shiatsu Instrukcja obsługi .....	51
<b>EN</b> Shiatsu foot bath Instructions for use.....	11	<b>NL</b> Shiatsu-voetenbad Gebruiksaanwijzing .....	58
<b>FR</b> Bain de pieds shiatsu Mode d'emploi.....	17	<b>DA</b> Shiatsu fodspa Betjeningsvejledning.....	65
<b>ES</b> Hidromasaje pies shiatsu Instrucciones de uso.....	24	<b>SV</b> Shiatsu-fotbad Bruksanvisning.....	71
<b>IT</b> Idromassaggiatore plantare shiatsu Istruzioni per l'uso .....	31	<b>NO</b> Shiatsu-fotbad Bruksanvisning.....	77
<b>TR</b> Shiatsu ayak banyosu Kullanım kılavuzu.....	38	<b>FI</b> Shiatsu-jalkakylpy Käyttöohje.....	83
<b>RU</b> Гидромассажная ванна для ног шиatsu Инструкция по применению .....	44		

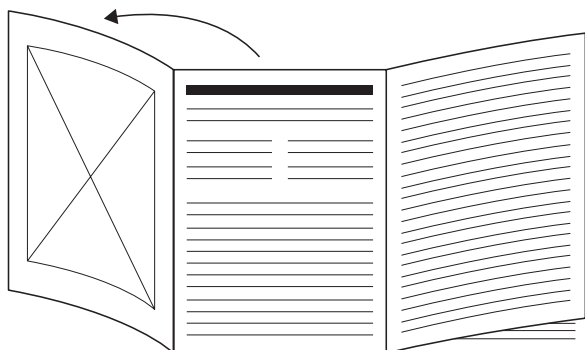- DE** Klappen Sie vor dem Lesen der Gebrauchsanweisung die Seite 3 aus.
- EN** Unfold page 3 before reading the instructions for use.
- FR** Dépliez la page 3 avant de lire le mode d'emploi.
- ES** Despliegue la página 3 antes de leer las instrucciones de uso.
- IT** Prima di leggere le istruzioni per l'uso aprire la pagina 3.
- TR** Kullanım kılavuzunu okumadan önce 3. sayfayı açın.
- RU** Перед чтением инструкции по применению разложите страницу 3.
- PL** Przed przeczytaniem instrukcji obsługi otworzyć stronę 3.
- NL** Vouw pagina 3 uit voordat u de gebruiksaanwijzing gaat lezen.
- DA** Fold side 3 ud, før du læser betjeningsvejledningen.
- SV** Vik ut sid. 3 innan du läser bruksanvisningen.
- NO** Åpne side 3 før du leser bruksanvisningen.
- FI** Käännä sivu 3 auki ennen käyttöohjeen lukemista.

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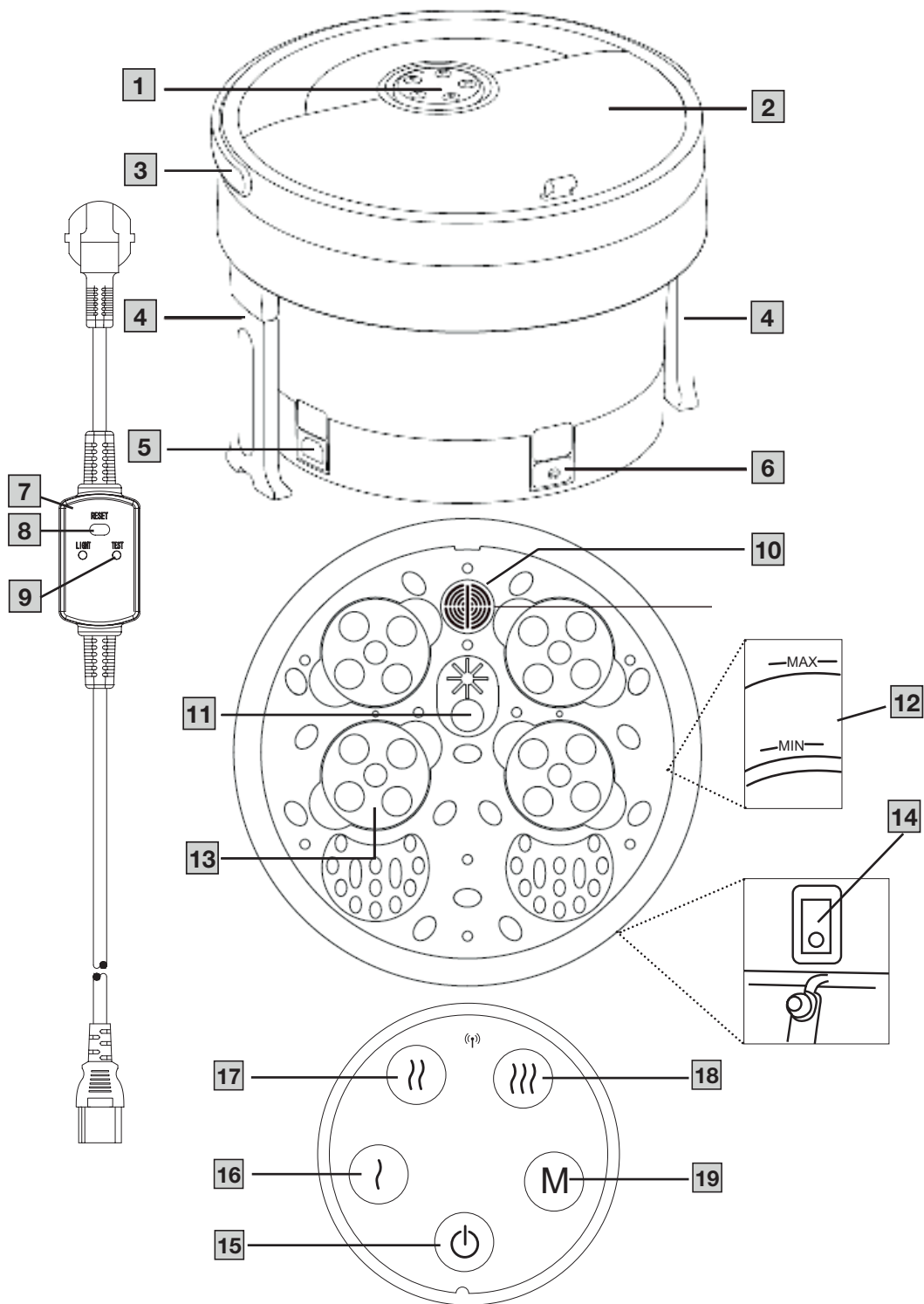**Read these instructions for use carefully. Follow the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.**

## ▲ WARNING

- The device is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 years and by people with reduced physical, sensory or mental capabilities or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the risks of using the device.
- Children must not play with the device.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- If any water leaks out of the device, it should no longer be used.
- The device has a hot surface. People who are not sensitive to heat must be careful when using the device.
- Under no circumstances should you open or repair the device yourself, otherwise proper functionality can no longer be guaranteed. Failure to comply with this instruction will void the guarantee.
- Cleaning and user maintenance must not be performed by children unless supervised.

## CONTENTS

1. Included in delivery .....	12
2. Signs and symbols .....	12
3. Warnings and safety notes .....	12
4. Device description.....	13
5. Usage .....	14
5.1 Please note the following before usage!.....	14
5.2 Preparing for usage.....	14
5.3 Starting usage.....	14
5.4 After usage.....	15
6. Cleaning and maintenance.....	15
6.1 Cleaning after each use .....	15
6.2 Thorough cleaning (weekly) .....	16
6.3 Storage.....	16
7. What if there are problems? .....	16
8. Accessories and replacement parts.....	16
9. Disposal.....	16
10. Technical specifications .....	16
11. Guarantee.....	16

## 1. INCLUDED IN DELIVERY


Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are included in the delivery. Before use, ensure that there is no visible damage to the device and accessories and/or replacement parts and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 Shiatsu foot bath
- 1 remote control
- 4 Shiatsu massage attachments
- 5 replacement filters
- 1 cover
- 2 stability supports
- 1 mains cable
- 1 copy of these instructions for use

## 2. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

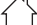	<b>WARNING</b> Indicates a potentially impending danger. If it is not avoided, there is a risk of death or serious injury.
	<b>CAUTION</b> Indicates a potentially impending danger. If it is not avoided, slight or minor injuries may result.
	<b>NOTICE</b> Indicates a potentially harmful situation. If it is not avoided, the system or something in its vicinity may be damaged.
	<b>Product information</b> Note on important information
	Read the instructions
	The electronic device must not be disposed of with household waste
	Manufacturer
	<b>CE labelling</b> This product satisfies the requirements of the applicable European and national directives.
	Marking to identify the packaging material. A = material abbreviation, B = material number: 1-7 = plastics, 20-22 = paper and cardboard
	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	Protection class II device
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU
	Protected against vertically falling drops of water
	Alternating current

	Direct current
	China Compulsory Certification
	For indoor use only
	United Kingdom Conformity Assessed Mark
	Importer

## 3. WARNINGS AND SAFETY NOTES

### Intended Purpose

#### WARNING

- This foot bath is intended to be used exclusively for the care of your feet.
- Never use the foot bath without water, as there is a risk of overheating.
- Do not use the foot bath with oily aromatic additives or oily bath salts.
- Do not use the foot bath if you suffer from a medical abnormality or injury in the area of your feet (e.g. open wounds, warts, athlete's foot).
- Do not use the foot bath while sleeping.
- Do not use the foot bath on animals.
- Never use the foot bath for more than 60 minutes (risk of overheating). Allow the foot bath to cool down for at least 15 minutes before using it again.
- Consult your doctor before using the foot bath if:
  - You are unsure whether the foot bath is suitable for you.
  - You suffer from a severe illness or have recently had surgery on your feet.
  - You suffer from diabetes or thrombosis, leg or foot conditions, e.g. varicose veins or venous inflammation, or pain of an unknown cause.
- The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

Note that the foot bath is suitable for people up to a shoe size of approx. EU 45.

### General safety notes

#### WARNING

- Danger of suffocation! Keep packaging material away from children.
- Do not pull, twist or bend the mains cable. Do not stick any needles or sharp objects into the device.
- Make sure that the mains cable is not routed or pulled over pointed or sharp objects.
- The device must not be connected to the mains during filling.
- Insert the mains plug into a socket that is protected against humidity, moisture and water splashes.
- Do not place the foot bath in a shower tray, bath tub or sink to fill it.
- If the device has been dropped or has suffered any other damage, it must no longer be used.
- Only specialist personnel may perform repairs on electrical devices. Improper repairs may put users in significant danger. In the event of a fault or damage, have the device repaired by a certified specialist workshop.
- Should a fault occur during use, switch the device off immediately and remove the mains plug from the socket.

- Always switch the device off and unplug it before moving or cleaning it.
- Check the foldable foot bath wall for possible cracks or water leaks before each use.
- Do not store sharp objects in the foot bath. Otherwise, this could damage the foot bath wall.
- Do not place feet or objects on the filled foot bath. This weight might cause the foot bath to fold up and water could escape.
- Danger of slipping! During use, water may escape from the foot bath over the rim.
- Do not step into the foot bath.
- Dispose of defective and discharged batteries immediately and properly (see chapter on disposal).
- Do not use modified or damaged batteries.
- Always select the correct battery type.
- Always insert the batteries correctly, taking into account the polarity (+/-).
- Never mix batteries of different manufacturers, capacities (new and used), size, or type within a device.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- If fluid from a battery comes into contact with your skin or eyes, wash the affected areas with water and seek medical assistance.
- **Choking hazard!** Keep batteries out of the reach of children. Seek medical attention immediately if swallowed. Swallowing them may cause burns, severe internal injuries, and death.
- Never allow children to replace batteries without adult supervision.

## Risk of electric shock

### ⚠ WARNING

- Make sure that the mains plug and mains cable do not come into contact with water or other liquids.
- Only use the foot bath in dry indoor areas. Never use the foot bath in the bath or a sauna.
- Only use the foot bath with dry hands.
- Never submerge the foot bath under water. Never pick up a foot bath that has been submerged in water. Disconnect the mains plug immediately.
- Only use the foot bath with the mains voltage specified on the device.
- Do not operate the foot bath during a storm.
- In the event of defects or malfunctions, switch the foot bath off immediately and disconnect it from the mains socket. Do not pull on the mains cable or on the foot bath to disconnect the mains plug from the socket.
- Never hold or carry the device by the mains cable. Keep the mains cable away from hot surfaces.
- Never connect the foot bath to the power cable if your feet are already in the water. Damage to the device can lead to fatal electric shocks.

## Risk of fire

### ⚠ WARNING

- Do not expose the device to high temperatures.
- Never leave the foot bath unattended when it is in operation, particularly if there are children nearby.
- Never operate the foot bath underneath a cover, such as a blanket, pillow, etc.
- Never use the foot bath near petrol or other highly flammable substances.
- Never use the foot bath on high-pile carpeting because the vents may be covered.

## Notes on handling batteries

### ⚠ WARNING

- **Risk of explosion! Risk of fire!** Failure to comply with the following points can result in personal injury or cause overheating, leakage, venting, breakage, explosion, or fire on the battery.
- This device contains non-rechargeable batteries which must not be charged.
- Do not throw batteries into a fire.
- Never charge, forcibly discharge, heat, disassemble, open, crush, deform, encapsulate, or modify batteries.
- Never short-circuit batteries or battery compartment contacts.
- Protect the batteries from direct sunlight, rain, heat, and water.
- Exposure of batteries to an environment with extremely high temperatures or an extremely low air pressure may result in explosion or leakage of flammable liquids and gases.

### ⚠ CAUTION

- Store batteries in a well-ventilated, dry, and cool place in a non-conductive container in which the batteries cannot be short-circuited to each other or by other metal objects.
- Keep batteries clean and dry.
- Keep batteries away from water.
- If the device is not going to be used for a long period of time, remove the batteries from the battery compartment.

### NOTICE

Do not use rechargeable batteries.

## 4. DEVICE DESCRIPTION

The associated drawings are shown on page 3.

- 1 Remote control
- 2 Cover
- 3 Carrying handle
- 4 Stability supports
- 5 Power connector
- 6 Substitute button
- 7 Mains adapter
- 8 Reset button
- 9 Test button
- 10 Filter
- 11 Indicator light
- 12 Fill-level mark
- 13 Detachable massage heads
- 14 Draining outlet
- 15 ON/OFF button
- 16 Heat 1 button
- 17 Heat 2 button
- 18 Heat 3 button
- 19 Massage button

## 5. USAGE

### 5.1 Please note the following before usage!


#### ⚠ WARNING

- **Never switch on the device without water in the basin.**
- **Never use the device with oily aroma additives or bath salts.**
- Never shave in the foot bath. The filter **10** could become clogged.
- Clean feet thoroughly before use.
- Clean the device after each use to remove hair from the filter **10**.
- If the remote control falls into water, remove the remote control and remove the battery. Allow the remote control to dry completely before using it again.


**i** The 30-minute automatic switch-off starts once a water temperature setting (1, 2 or 3) is selected by pressing a button and the corresponding water temperature has been reached. If you then switch to a different water temperature setting, the automatic switch-off will not be reset.

### 5.2 Preparing for usage


1. Remove the packaging.
2. Check the complete scope of delivery (including packaging) for damage.
3. Fold open the folded foot bath by pressing against the base of the foot bath with one hand and pulling the upper edge of the foot bath upwards with the other hand, for example.

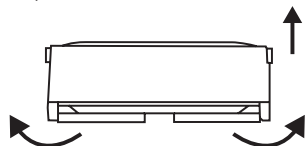Insert the two **4** stability supports into the side of the device until you hear and feel them click into place.

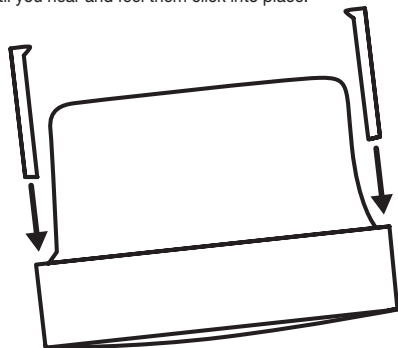4. Before plugging in the foot bath, always make sure that it is filled with water up to the minimum fill level mark "MIN". The water level must not exceed the maximum fill level mark "MAX" (10.5 cm); otherwise, water may overflow.
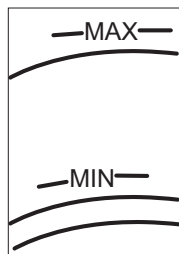
#### NOTICE

Make sure that the water temperature does not exceed 48°C, otherwise the seals may be damaged.

5. Place the foot bath on a firm and level surface.
6. Ensure that the mains cable does not pose a trip hazard.
7. Insert the mains adapter **7** into the power connector **5** on the foot bath.
8. Insert the mains plug into a suitable socket. The indicator light **11** starts to flash white. The foot bath is in standby mode.
9. Sit in a comfortable position. **Never use the foot bath standing up.**
10. Place your feet in the foot bath. Rubber feet on the base of the foot bath stop it from slipping.

**i** Do not use extension cords or multiple socket outlets as these may reduce the effectiveness of the protection and pose a safety risk.

Perform a functional test of the mains adapter at regular intervals. To do so, proceed as follows:

#### Resetting the mains adapter/function test

1. Plug the mains adapter **7** into a suitable socket.
2. Press the reset button **8** on the mains adapter. You should hear an audible "click". This indicates that the internal switch has been activated and power is applied. Observe the indicator light: If it flashes, a fault may have occurred. If it lights up continuously, the connection is normal.
3. Then press the test button **9**. You should hear another "click". This indicates that the power supply has been interrupted. This confirms that the residual current device is working properly.

#### Switching on the mains adapter

- Plug the properly reset mains adapter **7** firmly and directly into a fixed wall outlet.

### 5.3 Starting usage

**i** The 30-minute automatic switch-off starts once a water temperature setting (1, 2 or 3) is selected by pressing a button and the corresponding water temperature has been reached. If you then switch to a different water temperature setting, the automatic switch-off will not be reset.

You can operate the foot bath directly using either the remote control or the substitute button **6** on the foot bath. Both variants are presented in the following chapter.

## Usage with remote control

1. Make sure that the foot bath is correctly plugged in and that the indicator light **11** flashes white.
2. Press the ON/OFF button **15** on the remote control to switch the foot bath on. The indicator light **11** now pulsates white. The foot bath is now ready for use.
3. Press the massage button **19** on the remote control to switch on the massage. The massage will start (without heat). The indicator light **11** will continuously light up white. Press the massage button **19** on the remote control again to switch off the massage.
4. The foot bath also has three different heat settings (40°C, 43°C and 46°C). You can adjust the heat using the heat buttons **16** **17** **18** on the remote control.

	Heat 1 – button <b>16</b> Water temperature: 40°C	The indicator light <b>11</b> pulsates blue while the device is heating up. As soon as the water reaches 40°C, the indicator light <b>11</b> will light up blue continuously.
	Heat 2 – button <b>17</b> Water temperature: 43°C	The indicator light <b>11</b> pulsates purple while the device is heating up. As soon as the water reaches 43°C, the indicator light <b>11</b> will light up purple continuously.
	Heat 3 – button <b>18</b> Water temperature: 46°C	The indicator light <b>11</b> pulsates red while the device is heating up. As soon as the water reaches 46°C, the indicator light <b>11</b> will light up red continuously.

 You can combine the massage with heat. Alternatively, you can also activate the massage or heat individually.

 Different water temperatures can have a positive effect on how your body feels. Choose the right water temperature for your wellness experience:

- **Cold water** (24–29°C – without water heating) can have a refreshing and vitalising effect on the body and help you feel full of energy and vigour.
- **Lukewarm water** (30–40°C) supports an optimal wellness experience and can have a nurturing and reassuring effect.
- **Warm water** (approx. 42–48°C) can have a particularly soothing and warming effect on the body and, as a result, can support the immune system. Note that the water temperature should increase slowly and should not be set too high initially.

5. To switch off the foot bath before the 30-minute automatic switch-off time has elapsed, press the ON/OFF button **15** on the remote control. The foot bath switches to standby mode and the indicator light **11** starts to flash white. The foot bath switches off automatically after 20 minutes in standby mode.

## Usage with substitute button (button on device)

1. Make sure that the foot bath is correctly plugged in and that the indicator light **11** flashes white.
2. To switch the foot bath on, press and hold the substitute button **6** on the foot bath for 2 seconds. The indicator light **11** continuously lights up white and the massage starts (without heat).
3. The foot bath also has three different heat settings (40°C, 43°C and 46°C). Select the desired setting by pressing the substitute button **6** briefly.

1x short 	Water temperature: 40°C	The indicator light <b>11</b> pulsates blue while the device is heating up. As soon as the water reaches 40°C, the indicator light <b>11</b> will light up blue continuously. Massage is active.
2x short 	Water temperature: 43°C	The indicator light <b>11</b> pulsates purple while the device is heating up. As soon as the water reaches 43°C, the indicator light <b>11</b> will light up purple continuously. Massage is active.
3x short 	Water temperature: 46°C	The indicator light <b>11</b> pulsates red while the device is heating up. As soon as the water reaches 46°C, the indicator light <b>11</b> will light up red continuously. Massage is active.
4x short 	Water temperature: No heat	Starts over from the beginning. The indicator light <b>11</b> will light up white continuously. Massage is active without heat.

4. To switch the foot bath off before the 30-minute automatic switch-off time has elapsed, press and hold the substitute button **6** on the foot bath for 2 seconds. The foot bath switches to standby mode and the indicator light **11** starts to flash white.

## 5.4 After usage

1. Switch off the foot bath after every use. First, unplug the mains adapter **7** from the socket and then disconnect it from the power connector **5** of the foot bath. The indicator light **11** goes out.
2. Clean the foot bath as described in chapter “6. Cleaning and maintenance”.


## 6. CLEANING AND MAINTENANCE

### 6.1 Cleaning after each use

#### CAUTION

Make sure that no water gets inside the device! Make sure that the foot bath is switched off and unplugged every time before cleaning. Keep the mains plug and mains cable away from water!

1. Unplug the device and empty the water through the draining outlet after each use.
2. When unplugging, make sure that you first unplug the mains cable from the socket and then remove the cable from the device.
3. Empty all water from the tank (including any water remaining in the pipes).
4. Rinse the device with tap water.
5. To remove any remaining water, turn the foot bath over after draining the water from the draining outlet and place it upside down on a towel or sink to allow the remaining water to drain out.
6. Wipe the device with a soft cloth. Do not use steel wool.

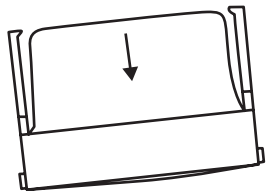
## 6.2 Thorough cleaning (weekly)

### CAUTION

Make sure that no water gets inside the device! Make sure that the foot bath is switched off and unplugged every time before cleaning. Keep the mains plug and mains cable away from water!

- Remove the massage heads **13** from the device for thorough cleaning. To remove the massage heads **13**, first press the middle button (A) and then pull the massage head upwards.

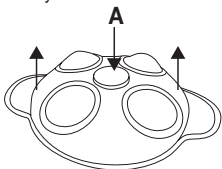- Wash the device with a mild household cleaner (not foaming, e.g. neutral vinegar).
- Remove the filter **10**. Wash the filter carefully by hand with warm water and a neutral cleaning product (do not use brushes).
- Allow the filter **10** to air dry completely. The filter should be checked regularly and replaced after 12 months. Replace immediately in case of severe yellowing, hardening or damage. 5 replacement filters are included in the package.

## 6.3 Storage

- If you do not intend to use the device for an extended period of time and it is completely dry, fold it up.
- Cover the device with the cover **2**.
- Place the device in its original packaging and store it in a safe, well-ventilated, dry and cool place.

## 7. WHAT IF THERE ARE PROBLEMS?

Problem	Possible cause	Solution
Unusual noises	Water level below the "MIN mark"	Check the water level and fill up to the MIN mark.
Water no longer flows	Nozzles are clogged	Disconnect the device from the mains and carefully remove any visible blockages from the nozzles using a toothpick or a soft brush.
Water no longer flows	Filter dirty/full	Remove the filter and insert a new one. 5 replacement filters are included in the package.
The device emits 5 short acoustic signals. None of the buttons respond any more.	Temperature sensor faulty	Remove your feet from the foot bath. All functions are deactivated. Switch the device off and on again by unplugging it.
Device emits 10 short acoustic signals, indicator light flashes red 10 times. None of the buttons respond any more.	Water temperature exceeds 50°C	Remove your feet from the device, remove the water and switch the device off and on again by unplugging it.

## 8. ACCESSORIES AND REPLACEMENT PARTS

Article	Article number
Support feet	110314
FB 55 adapter	110315
FB 55 UK adapter	110317
FB 55 remote control	110316

## 9. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Observe the local regulations for material disposal. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal. You can obtain the location of collection points for old devices from the local or municipal authorities, local waste disposal companies or your retailer, for example.

## 10. TECHNICAL SPECIFICATIONS

Mains operation:	Input 220-240 V, 50/60 Hz, 800 W Output 3.0 V $\approx$ 12 mA
Size:	$\varnothing$ approx. 380 mm * 280 mm
Weight (without cable):	approx. 4.1 kg
Power consumption of device in off mode:	0.0 W
Power consumption of device in standby mode:	0.5 W
Period after device reaches automatically standby mode:	30 minutes

## 11. GUARANTEE

Further information on the guarantee and guarantee conditions can be found in the guarantee leaflet supplied.

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BEURER UK Ltd.
Beurer GmbH • Söflinger Str. 218 • 89077 Ulm, Germany