# chal-tec gmbh

elektronik  star

## Tageslichtlampe

10006251

elektronik  star  
[www.elektronik-star.de](http://www.elektronik-star.de)

Telefon: 01805-019 180\* (Mo-Fr 10-13 und 14-18 Uhr)

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## Normative references

<b>Guidance and manufacturer's declaration – electromagnetic immunity</b>			
The "Daylight" is intended for use in the electromagnetic environment specified below. The customer or the user of the "Daylight" should ensure that it is used in such an environment.			
<b>Immunity test</b>	<b>IEC 60601 test level</b>	<b>Compliance level</b>	<b>Electromagnetic environment – guidance</b>
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % $U_T$ (>95 % dip in $U_T$ ) for 0,5 cycle 40 % $U_T$ (60 % dip in $U_T$ ) for 5 cycles 70 % $U_T$ (30 % dip in $U_T$ ) for 25 cycles <5 % $U_T$ (>95 % dip in $U_T$ ) for 5 sec	<5 % $U_T$ (>95 % dip in $U_T$ ) for 0,5 cycle 40 % $U_T$ (60 % dip in $U_T$ ) for 5 cycles 70 % $U_T$ (30 % dip in $U_T$ ) for 25 cycles <5 % $U_T$ (>95 % dip in $U_T$ ) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the "Daylight" requires continued operation during power mains interruptions, it is recommended that the "Daylight" be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE $U_T$ is the a.c. mains voltage prior to application of the test level.			

## Guidance and manufacturer's declaration – electromagnetic immunity

The "Daylight" is intended for use in the electromagnetic environment specified below. The customer or the user of the "Daylight" should ensure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 V	Portable and mobile RF communications equipment should be used no closer to any part of the "Daylight", including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.  <b>Recommended separation distance</b>  $d=1.2\sqrt{P}$ $d=1.2\sqrt{P}$ 80MHz to 800MHz $d=2.3\sqrt{P}$ 800MHz to 2.5 GHz  where $P$ is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and $d$ is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, <sup>a</sup> should be less than the compliance level in each frequency range. <sup>b</sup> Interference may occur in the vicinity of equipment marked with the following symbol:  
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

<sup>a</sup> Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the "Daylight" is used exceeds the applicable RF compliance level above, the Medical Daylight should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the "Daylight".

<sup>b</sup> Over the frequency range 150 kHz to 80 MHz, field strengths should be less than  $[V_1]$  V/m.

<b>Guidance and manufacturer's declaration – electromagnetic emissions</b>		
The "Daylight" is intended for use in the electromagnetic environment specified below. The customer or the user of the "Daylight" should ensure that it is used in such an environment.		
<b>Emissions test</b>	<b>Compliance</b>	<b>Electromagnetic environment – guidance</b>
RF emissions CISPR 11	Group 1	The "Daylight" uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The "Daylight" is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

<b>Recommended separation distances between portable and mobile RF communications equipment and the Medical Daylight</b>			
The "Daylight" is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Medical Daylight can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the "Daylight" as recommended below, according to the maximum output power of the communications equipment.			
<b>Rated maximum output power of transmitter W</b>	<b>Separation distance according to frequency of transmitter m</b>		
	<b>150 kHz to 80 MHz</b> $d = \left[ \frac{3,5}{V_1} \right] \sqrt{P}$	<b>80 MHz to 800 MHz</b> $d = \left[ \frac{3,5}{E_1} \right] \sqrt{P}$	<b>800 MHz to 2,5 GHz</b> $d = \left[ \frac{7}{E_1} \right] \sqrt{P}$
0,01	<b>0.12</b>	<b>0.12</b>	<b>0.23</b>
0,1	<b>0.38</b>	<b>0.38</b>	<b>0.73</b>
1	<b>1.2</b>	<b>1.2</b>	<b>2.3</b>
10	<b>3.8</b>	<b>3.8</b>	<b>7.3</b>
100	<b>12</b>	<b>12</b>	<b>23</b>
For transmitters rated at a maximum output power not listed above, the recommended separation distance $d$ in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where $P$ is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

## Safety Notes

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

### Warning

- Do not use in humid surroundings.
- Before use, ensure that all packaging materials are removed and that there is no visible damage to the unit or accessories. When in doubt, do not use the unit and contact your dealer or the customer service address provided.
- The "Daylight" is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the "Daylight".
- Keep the cord away from hot objects and open flames.
- The mandatory minimum distance between the "Daylight" and flammable objects is 1.0 meter.
- Connect the unit only to the mains voltage listed on the type plate.
- Always remove the mains plug and let the unit cool off before you touch it.
- Do not touch the unit with wet hands while it is plugged in; do not allow any water to be sprayed onto the unit. The unit must be operated only when it is completely dry.
- Protect the unit from strong impacts.
- Avoid pulling the cord when removing the plug from the socket. Instead, hold the plug.
- If the cord or the housing is damaged, contact customer service or the dealer, since special tools are needed to make repairs. Never try to repair the appliance yourself.
- Disconnection from the mains supply is only guaranteed if the mains plug has been removed from the socket; do not position the daylight so that it is difficult to operate the plug device.
- The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- The device requires no calibration.
- The device is not repairable and contains no user serviceable parts.
- The operator shall not touch Daylight and the patient simultaneously.
- No modification of this equipment is allowed.

## Notes

### General

#### Important

- This unit is not intended for commercial or clinical use, but only for private household use.
- If you are taking medication such as pain relievers, antihypertensives or antidepressants, consult your doctor before using the "Daylight".
- People with retinal diseases as well as diabetics need to be examined by an ophthalmologist before beginning to use the "Daylight".
- Please do not use if you have an eye disease such as cataracts, glaucoma, diseases of the optic nerve in general and inflammation of the vitreous humour.
- If you have health concerns of any kind, contact your general practitioner!

### Repairs

#### Important

- Only open the unit to replace the tube. Please do not attempt to repair the unit yourself. This could result in serious injuries. Failure to observe this regulation shall void the warranty.
- If this appliance is damaged in any way or it is not working, it should not be operated. Please try to fix the problem according to FAQ table. Never try to repair the appliance yourself.

## Unpacking and inspecting

Every attempt is made to ensure your order is accurate and complete. However, to be sure that your order is correct, verify the contents of the box against your packing slip.


The "Daylight" is designed for simplicity of operation and set-up and requires minimal assembly. The following items are included your box:

- 1 (one) "Daylight"
- 1 (one) Operator's Manual

Carefully inspect each item as it is unpacked for any signs of damage which may have occurred during shipment.

- Check the components according to the packing list.
- Check for any damage or defects. Do not attempt to setup the "Daylight" if anything is damaged or defective. Contact your supplier immediately if anything is found to be damaged or defective.

## Unit description

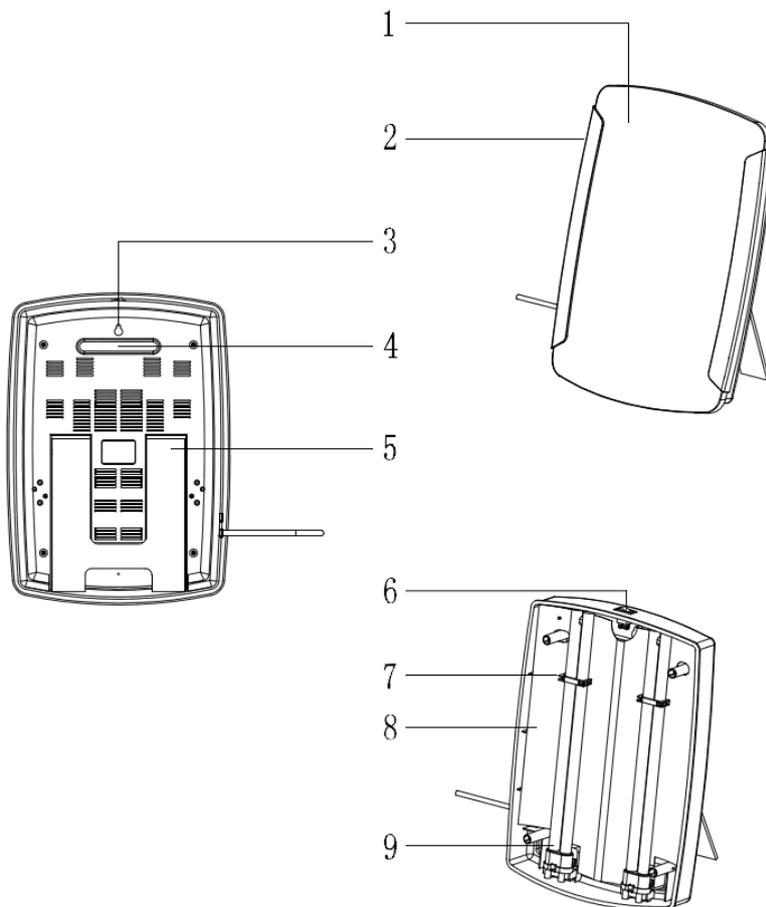1. **Screen**
2. **Brackets**
3. **Mounting hole**
4. **Carrying handle**
5. **Fold-out stand**
6. **On/off switch**
7. **Retaining clips**
8. **Reflector**
9. **Tubes, 2 x 36W**

## Preparing the appliance for use

### Setup

1. Remove the unit from the packaging.
2. Place the unit on a solid surface, such as a table, and adjust the folding tripod.
3. It is recommended that you keep your face 20 - 50 cm away from the light during use.

### Mains connection

Connect the unit only to the mains voltage listed on the type plate.

- Push the plug completely into the mains socket.

#### Note

- Ensure that there is a socket near the place where the lamp will stand.
- Lay the mains cable in such a way that no one can trip over it.

### Fastening of Daylight on the wall

There is a mounting hole on the back side of the Daylight. At this mounting hole the unit can be fixed on the wall with screw. Depending on the type of wall, an appropriate screw must be used, may be with pegs. E.g. for concrete wall, M4, 25mm long Stainless Steel Self-Tapping Screw which is screwed into a anchor. The screw and the anchor are not included in the scope of delivery.

#### Important

- If the mounting hole or back enclosure is cracked or damaged. Do not fix the Daylight on the wall. Please contact customer service or the dealer for assistance.

## Operation

### • Switching on the lamp

Press the ON/OFF button on the upper edge of the lamp to switch on the "Daylight".

#### Note

The brightness of the tubes may be visibly different when switched on. The tubes develop their full luminosity after a couple of minutes.

### • Taking the light

Now sit down in front of the light or as close to it as comfortable, between 20 cm and 50 cm. You can go about your normal activities while using it. You can read, write, make telephone calls, etc. When using the lamp, try to maintain the recommended distance of 20 - 50 cm between your face and the lamp.

Usage guidelines are as follows:

- Approx. 2 hours at a distance of 60cm
- Approx. 1 hour at distance of 40 cm
- Approx. 30 minutes at a distance of 20cm

In effect, the closer you are to the source of light, the less time you need to use it.

#### Note

- Do not look directly into the light continuously! Some eye pain or mild headaches may be experienced during initial sessions but these should pass. If in doubt just lower your treatment times until comfortable again.
- Do not leave the appliance unattended in the duration.
- Do not touch the device in the duration. Let the device cool down for at least 10 minutes before repositioning it.

- **Enjoying light over longer periods**

During the darkest part of the year, repeat the treatment for at least 7 consecutive days, or for longer, depending on your individual requirements.

- **Switching off the "Daylight"**

Press the ON/OFF button on the upper edge of the lamp to switch off the "Daylight". Pull the mains plug from the socket.

**Note**

- Let the appliance cool down first before storing it away.

## Replacing the lamp tubes

### Replacing the lamp:

In the event that the lights are faulty, or when they come to end of their life, the two light tubes can be replaced. NOTE: Do not replace the lamp tubes which are not provided by the supplier. To do so, proceed as follows:

**Warning**


- Always unplug the "Daylight" and let it cool off.
- For safety reasons, do not remove any components or screws other than those specified below.

1. Slacken the retaining screws on the back using a screwdriver, lay the apparatus on its back and remove the screen.

**Note**

The two brackets on the front of the screen will become loose when the screws are slackened. Take care not to lose them, or the screen cannot be refitted correctly.


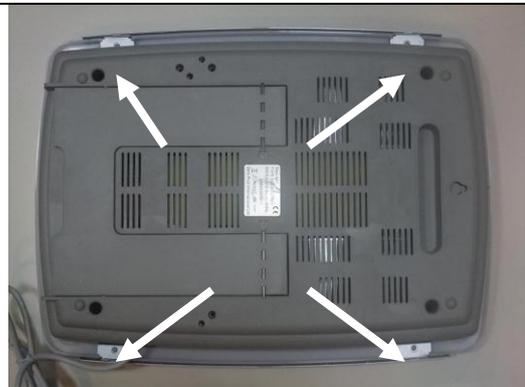
2. Bend the tube's four transparent retaining clips by pulling them gently upwards.

**Note**

Cable ties are used for transport protection only.
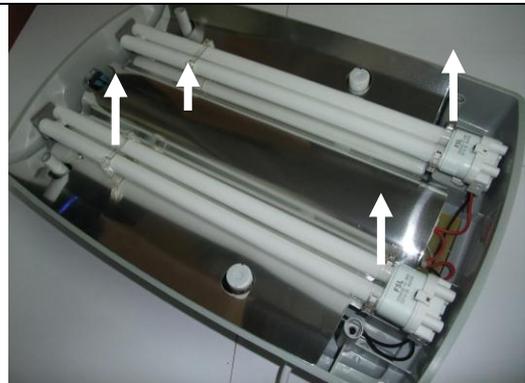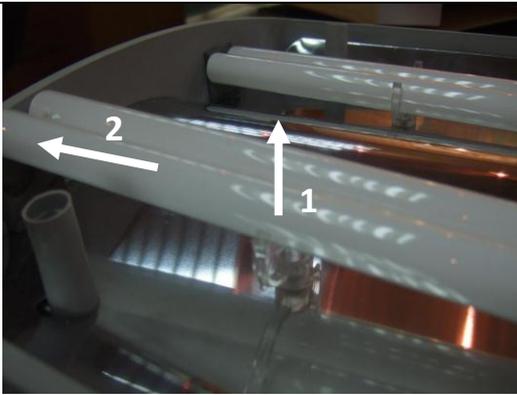
<p>3. Carefully pull the tube up over the housing rim (1) and out of the housing, then sideways out of the frame.</p>	
<p>4. When inserting the tube, following these instructions in reverse order.</p>	

### Cleaning Up a Broken Tube

The light tubes contain a small amount of mercury sealed within the glass tubing, when a tube breaks in your home, some of this mercury is released as mercury vapor. To minimize exposure to mercury vapor, we recommend that users follow the cleanup and disposal steps as follows:

#### Before Cleanup:

- have people and pets leave the room and air out the room for 5-10 minutes by opening a window or door to the outdoor environment.
- Shut off the central forced air heating/ air-conditioning system, if you have one.
- Collect material needed to clean up broken tube, just like: stiff paper or cardboard, sticky tape, damp paper towels or disposable wet wipes.

#### During Cleanup:

- Vacuuming is not recommended unless broken glass remains after all other cleanup steps have been taken. Vacuuming could spread mercury-containing powder or mercury vapor.
- Be thorough in collecting broken glass and visible powder.
- Place cleanup materials in a sealable container.

#### After Cleanup:

- Promptly place all bulb debris and cleanup materials outdoors in a trash container or protected area until materials can be disposed of properly. Avoid leaving any bulb fragments or cleanup materials indoors.
- If practical, continue to air out the room where the bulb was broken and leave the heating/air conditioning system shut off for several hours.

### Cleaning

#### **Important**

- Before you clean the unit, always switch it off, unplug it and let it cool off.
- Do not touch the unit with wet hands.
- Clean using a damp cloth, to which you can apply a little detergent if necessary.
- Do not let water run into the appliance or spill water onto the appliance.

## Storage

When the "Daylight" is not being used, it should be stored in a dry location out of reach of children. Do not keep the "Daylight" connected to a power source. The "Daylight" is designed to be stored in the following environmental conditions:

- 1) Ambient temperature from 41 Degrees Fahrenheit (5 centigrade) to 122 Degrees Fahrenheit (50 degrees centigrade).
- 2) Relative humidity of max. 90%.

## FAQ

Failure	Possible causes	Solutions
No light	Main Power OFF	Turn on the switch
	No Power	The appliance is not correctly connected to your power supply is defective.
	Loose Light tube	Re-install the light tube.
One side of the product is out of order	Defective light tube	Change the light tube
	Light tube at end of life	Change the light tube
	Loose Light tube	Re-install the light tube.

### Important

This chapter lists the questions most frequently asked about the appliance. If you cannot find the answer to your question, please contact customer service or the dealer for assistance..

## Technical specification

<b>Dimension (WHD)</b>	340 x 480 x 80 mm
<b>Weight</b>	2.0kg
<b>Lighting elements</b>	2 x 36W tubes (2G 11), Daylight (full spectrum)
<b>Nominal data for replacement lighting elements (provided by light tube supplier)</b>	PLL 36W 2G11 4pin H Philips - MASTER PL-L 36W/865/4P
<b>Fuse</b>	T3.15AH250V
<b>Power</b>	72 watts
<b>Light intensity</b>	> 10,000 lux in 15cm
<b>Mains connection</b>	220 V ~240V, 50 Hz
<b>Operating Environment</b>	50 – 77 °F/ 10 – 25 °C, and relative humidity from 30% to 90%. Pressure altitude: <2000m
<b>Storage Environment</b>	41 - 122°F/ 5 – 50 °C, and relative humidity of max. 90%. Pressure altitude: <2000m
<b>Transport Environment</b>	41 - 122°F/ 5 – 50 °C, and relative humidity of max. 90%. Pressure altitude: <2000m
<b>Product Classification</b>	<ol style="list-style-type: none"> <li>1. Class II equipment</li> <li>2. No applied part</li> <li>3. Protection against ingress of water: IPX0</li> <li>4. Not category AP / APG equipment</li> <li>5. Mode of operation: Continuous operation.</li> </ol>

Subject to technical modification.