## Lichtdusche

10026416

# KLARSTEIN

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Dear Customer,

Congratulations on purchasing this equipment. Please read this manual carefully and take care of the following hints to avoid damages. Any failure caused by ignoring the mentioned items and cautions mentioned in the instruction manual are not covered by our warranty and any liability.

## Technical Data

Item number	10026416
Power supply	220-240 V ~ 50-60 Hz
Light intensity (MAX)	10.000 Lux

## Safety Instructions

To reduce the risk of fire, electric shock and damage:

- Do not expose this device to rain.
- Do not place objects on the unit filled with liquids, such as vases.
- Use recommended accessories only.
- Do not repair this device yourself. Refer repairs to qualified service personnel.
- Do not drop any metal objects into this unit.
- Do not place heavy objects on top of this unit.

### Power cord / power supply

- Make sure that the power supply voltage corresponds to the specifications on this unit.
- Plug the power cord into the outlet securely.
- Do not pull the cable, do not bend it and do not place heavy objects on it.
- Do not handle the plug with wet hands.
- Grasp the plug when pulling out of the electrical outlet.
- Do not use a damaged power plug or a damaged electrical outlet.
- Install this unit so that the power cord can be unplugged immediately when a problem occurs.

### Small objects / Packaging Parts

Store small objects (eg screws and other mounting hardware, memory cards) and packaging parts out of the reach of children. Do not let you play with the film. Risk of suffocation!

### Location

- Set the appliance on a dry, flat, waterproof and heat-resistant surface.
- The socket must be readily accessible to disconnect the unit quickly in an emergency from the electrical circuit.

- Do not install or position this unit in a bookcase, built-in wardrobe or any other confined space. Ensure good ventilation of the unit.
- Do not place near to an amplifier or other equipment that may be hot. The heat can damage the unit.
- Do not expose this unit to direct sunlight, high temperatures, excessive moisture and excessive vibration.

### Power Supply

- The power adapter can operate on AC 220-240 V only. Use a suitable plug adapter if the plug does not fit into the power outlet.

### Intended Use

- This device is used for producing light. It is exclusively meant for this purpose and must only be used for this purpose. It may only be used in the manner as described in this manual.
- This product is not designed for use by persons (including children) who have reduced physical, sensory or mental capacities or persons with no experience or knowledge, unless they are instructed and supervised in the use of the appliance by a person responsible for their safety.

### Transport

- Please keep the original packaging. To an adequate protection during transport of the device, pack the unit in its original packaging.

### Cleaning

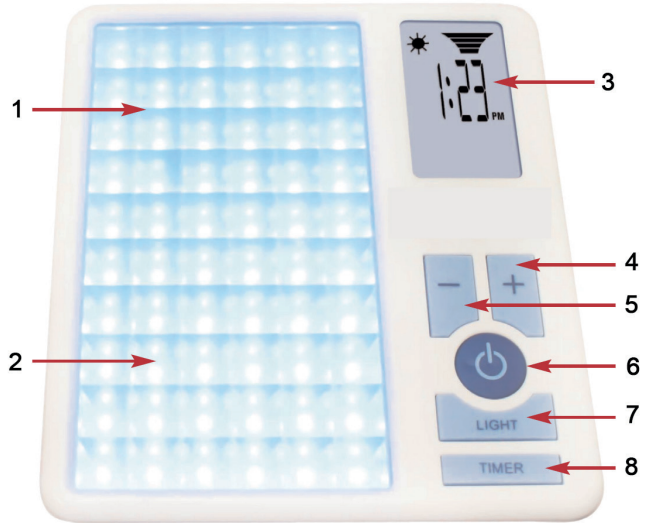
Do not use volatile liquids, such as insecticide spray. Don't clean the surface with strong pressure otherwise it could get damaged. Rubber or plastic should not be in contact with the device over a longer period. Clean the appliance and the floor using a soft cloth and warm water, if necessary, a soap solution. Never submerge the base in water or any other liquid.

### Repairs

- Repairs to the appliance should only be carried out by an authorized service center.
- A remodeling or modification of the affected product safety. Risk of injury!
- Open the product never arbitrarily and perform repairs never yourself!
- Incorrect reassembly may result in malfunction or total failure.
- Open in no case the unit. We assume no liability for consequential damages. The warranty expires at foreign interference with the device.

# Product Description

- 1 60 blue LEDs
- 2 Anti-glare-Reflektor/  
Diffuser
- 3 LCD Screen
- 4 [+] to increase the LED  
intensity and set the Timer
- 5 [-] to decrease the LED  
intensity and set the Timer
- 6 ON/OFF
- 7 Light
- 8 Timer: to set the  
treatment duration

## Display Indication

	Treatment session timer indicator
	Battery capacity / charge indicator
	LED intensity indicator
	Sun Symbol will always remain on

# Use and Operation

## Charge Battery

1. Connect the dc plug of adapter to the unit. Plug in the adapter to the wall socket.
2. The initial charge should be a minimum of 12 hours. Subsequent re-charging will take 5-6 hours.
3. Battery indicator on LCD screen will flash whilst charging and stop flashing indicating 3 bars when fully charged.

## Position your Klarstein Blue Sun Correctly

- You should sit in a comfortable chair in front of a table or desk.
- Place Blue Sun at an angle of around 45° to the left or right of your vision directly in front of you. (See diagram 1).
- Adjust stand so LED's are directed to your eye level.
- Blue Sun should be positioned NOT CLOSER than 50cm (20 inches), AND NOT FURTHER than 70cm (30 inches) from your eyes.
- The environment where you are using Blue Sun should be a well lit room – Natural daylight or electric lighting if it is dark.
- You can carry on doing normal activities whilst using Blue Sun ie. Reading, writing, using a computer or watching TV.

## Easy Operation

1. Press the round power button to turn on. LCD display will indicate 12:00.
2. To set the time press timer button using finger of one hand and hold down, whilst at the same time press + or – button with finger of other hand to adjust time.

If you keep + or – button depressed, display will move very quickly! You will not need to set time again even if unit is switched off. It will memorise the time set.

## Manual Operation

Press light button and Blue LED's will come on. Press + or – button to required level of intensity. There are 4 setting which are indicated on the LCD display. Press light button to switch off LED's after you have ended your session. Switch power button OFF and LCD display will disappear. (If you forget to turn power OFF, it will Auto switch off.)

## Timed Session

- Press timer button and timer symbol will appear on LCD screen with default time of 10 mins underneath which will flash.
- To adjust required length of session press + or – button. (Up to 60 mins) Blue LED's will automatically come on after 5 seconds.
- If you want to re-set session time again press light button so LED's go off. Press timer button again and the + or – button to select required session length. Session length time will count down on the LCD screen.
- LED's will go off after your selected session length has ended.
- Auto Memory The Blue Sun will automatically remember time, LED intensity and session length from your last session.

## Re-Set Button

The re-set button is located at the back of the unit. Insert tip of ball point pen and press until you here “click”. Unit will re-set to factory settings.

## When you should use Klarstein Blue Sun

- If you are feeling generally rather “down” and “lacking in energy” particularly during the dark winter months – you may benefit greatly from a daily session of around 30 minutes using SAD Solutions Blue Sun at 50% intensity. This will help compensate for the lack of natural sunshine.
- Light affects your circadian rhythm or body clock. In simple terms this means the body “wakes up and all functions are fully on” during daylight hours. During darkness the body “slows down” and all functions reduce until it goes into “sleeping mode”.
- This is why you may find it easier to wake up and start your day in summer months when it is light. Conversely in winter when it is dark, you may feel less inclined to wake up.
- With the stresses of modern day life, your natural circadian rhythm and body clock can go out of sync.
- It is possible to use Klarstein Blue Sun to re-set your body clock to suit your routine and lifestyle. Below are some examples of how Klarstein Blue Sun can be very effective.
- It is recommended that you start using Blue Sun at 50% intensity and gradually increase to your comfort level.

## Waking Up and Sleeping

- If you fall asleep in front of the TV and when you go to bed have difficulty sleeping – then use Blue Sun for 30 mins while you are watching TV and you are less likely to go to sleep in your armchair and more likely to sleep when you go to bed.
- If you find yourself going to bed early, and then waking up too early in the morning, then also use Blue Sun to keep you awake longer in the evening. You are then more likely to wake up later.
- If you find it difficult to wake up in the morning, and feel “slow and drowsy” then try a 30 minute session of Blue Sun when you wake up – then you will feel refreshed and ready to start your day.

## Shift Work

- If you have trouble staying awake particularly at the start of a night shift routine then use Blue Sun for 30 minutes before you start work and repeat for another 15 minutes during the night if you start to feel drowsy.
- This should help your body clock adjust to this new routine. If possible avoid bright light when coming off a night shift and preferably wear sun glasses. You may find this helps you go to sleep when you return home.
- After coming off a night shift schedule, the temptation is to fall asleep because your body clock is used to this routine. Use Blue Sun to keep you awake until evening – and again this will help re-set your body clock so you sleep at night rather than during the day.

## Jet Lag

Long distance air travel East or West across different time zones can play havoc with your body clock. As a general tip try planning your arrival time at your destination in the morning rather than the evening. Sleep on the plane if you are able. Then try to keep awake for as long as possible after you reach your destination, preferably in natural daylight – and using your Blue Sun can be a great help in keeping you awake until it gets dark.

## Cleaning

Use a soft dry cleaning cloth only. Do not use a damp cloth, or any cleaning solution. Take care not to scratch the lens.

## Hints on Disposal

According to the European waste regulation 2002/96/EG this symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. For more detailed information about recycling of this product, please contact your local council or your household waste disposal service.

## Declaration of Conformity

Producer: Chal-Tec GmbH, Wallstraße 16, 10179 Berlin, Germany.

This product is conform to the following European Directives:  
2014/30/EU (EMC)  
2014/35/EU (LVD)  
2011/65/EU (RoHS)

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Fabricant : Chal-Tec GmbH, Wallstraße 16, 10179 Berlin, Allemagne.