Niwe

10030232 10030233 10030234

KLAR  FIT

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2004/108/EC (EMC)  
2011/65/EU (RoHS)

Dear Customer,

Congratulations on purchasing this product. Please read and follow these instructions, in order to avoid damaging the item. We do not cover any damages that may arise from improper use of the item or the disregard of the safety instructions.

## General Safety Notice

- Read all instructions before using.
- Do not place on or near a hot gas or electric burner or in a heated oven.
- Do not expose the appliance to extreme temperatures. The appliance shall be used in moderate climates.
- Do not expose to water or moisture. Clean with a dry cloth.
- The use of attachments that are not recommended by the appliance manufacturer may result in personal injury.
- The appliance is not a toy. Do not let the children play with it. Never let children insert foreign objects into the appliance.
- Do not use the appliance for other than intended use. This appliance has been designed solely for domestic use; the manufacturer's liability shall not be engaged in the event of any professional use, by any misuse or by any use not complying with the operating instructions.
- Save these instructions.
- This machine is not suitable for therapeutic purposes.
- People such as children and handicapped persons should only use the machine in the presence of another person who can give aid and advice.
  
- Never leave the device unattended when plugged in. Unplug from the outlet when not in use and before wiping.
- Never operate this appliance if cord or plug are damaged, if it is not working properly or has been dropped into water.
- Never operate the appliance with the air openings blocked. Keep any openings free of lint, hair etc.
- Do not use outdoors.
- Never drop or insert any object into any opening.
- Do not operate where aerosol spray products are being used or where oxygen is being administered.
- Do not stand on the appliance. Use only while seated.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation.
  
- Remove all scarves, neckties, necklaces and jewelry before using this massager. Any individual who may be pregnant, has a pacemaker, suffers from diabetes, phlebitis and/or thrombosis, is at an increased risk of developing blood clots, or who has pins, screws, artificial joints or other medical devices implanted should consult a physician before using this device.
  
- Do not use on an infant or invalid or on a sleeping or unconscious person.
  
- Never remove the cover.
  
- **Small objects:** Keep small objects (i.e. screws, mounting material, memory cards etc.) and packaging out of the range of children. Do not let Children play with foil. Choking hazard!
  
- **Transporting the device:** Please keep the original packaging. To gain sufficient protection of the device while shipment or transportation, make sure to put it back in its original packaging.
  
- **Cleaning the surface:** Do not use any volatile liquids, detergents or the like. Use a clean dry cloth.

# Overview

adapter


Office Chair

Sofa

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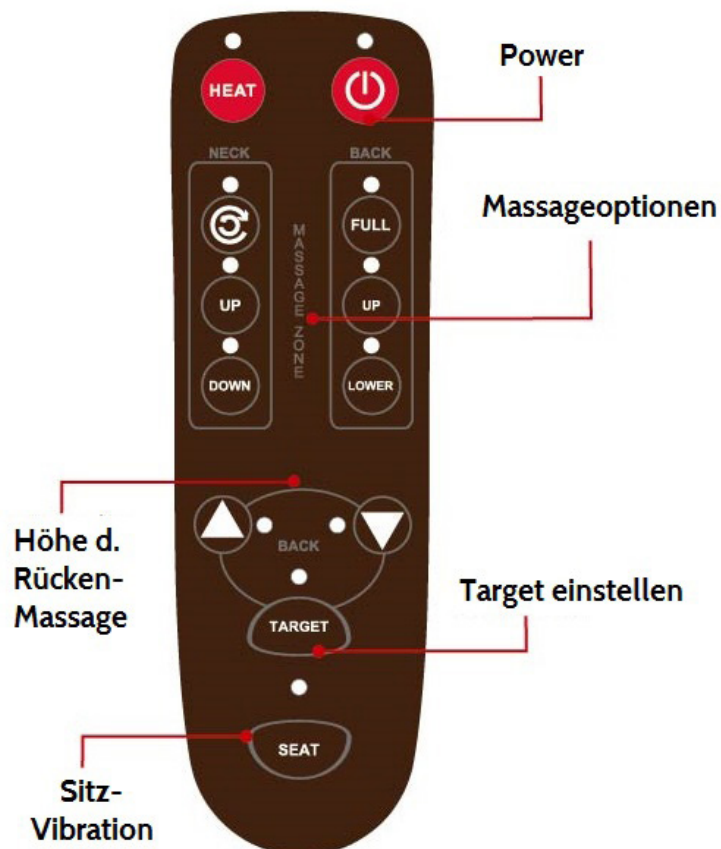## Operation

**Preparation:** Place the massager on a sturdy and comfy chair. Fix the back strap to your chair/seat. Plug the AC adapter into a power supply and insert the DC end into the dual-use plug of the cushion.

**Operation:** Use the RC for all functions (see above). You can select and combine the massage programs according to your preference and during operation. After the massage is completed, the rollers will return to storage position. Unplug from the electric outlet when unit is not in use. Please note: The unit shuts off automatically after approx. 15 minutes when not used.

**Troubleshooting:** Massager does not start? Ensure that power is on, plug is connected and voltage is adequate. After prolonged use, the motors may need to cool down to protect them. Please wait 30 minutes before trying again. Massage rollers stop during massage session? The mechanism may be stalled when subjected to excessive force. Frictional sound? May be a result of the backrest cover. Motor and air sounds may occur.

## Environment Concerns

According to the European waste regulation 2002/96/EG this symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. For more detailed information about recycling of this product, please contact your local council or your household waste disposal service.

## Declaration of Conformity

**Producer: CHAL-TEC GmbH, Wallstraße 16, 10179 Berlin, Germany.**

This product is conform to the following European directives:

2004/108/EC (EMC)  
2011/65/EU (RoHS)