Massagepistoie

Massage gun

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Pistolet de massage

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10038969

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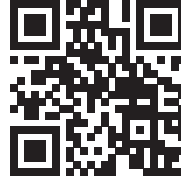
**Dear customer,**

Congratulations on the purchase of your appliance. Please read the following instructions carefully and follow them to prevent potential damage. We accept no liability for damage caused by disregarding the instructions or improper use. Please scan the QR code to access the latest operating instructions and for further information about the product.

---

**CONTENTS**

---

Safety Instructions	10
Technical data	11
Scope of delivery	11
Essays	12
Control panel and displays	13
Operation	14
Disposal Considerations	14

---

**MANUFACTURER & IMPORTER (UK)**

---

**Manufacturer:**

Chal-Tec GmbH, Wallstraße 16, 10179 Berlin, Germany.

**Importer for Great Britain:**

Chal-Tec UK limited  
Unit 6 Riverside Business Centre  
Brighton Road  
Shoreham-by-Sea  
BN43 6RE  
United Kingdom

---

## SAFETY INSTRUCTIONS

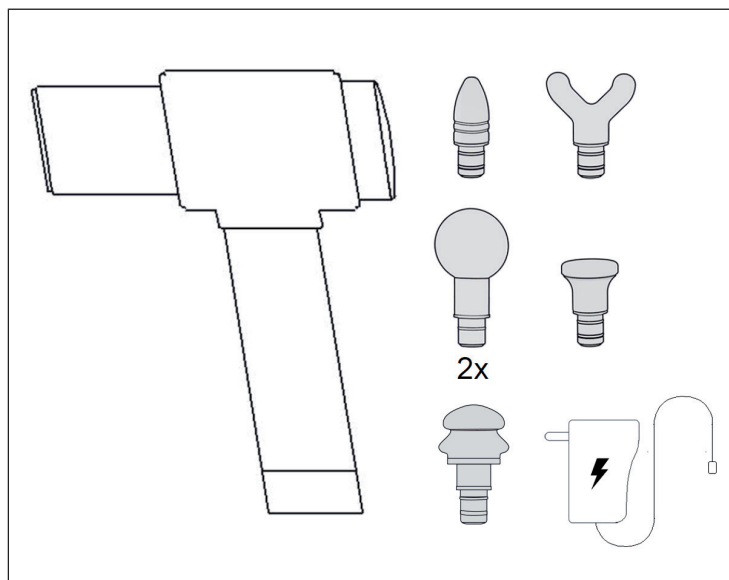
---

- This appliance may only be used by adults.
- Discuss the use of this device with your doctor, especially if the following applies to you:
  - You are an osteoporosis patient.
  - You have heart problems.
  - You wear an electronic medical device.
  - You are pregnant.
  - You are at risk of thrombosis or have thromboses.
  - You have a body trauma or other medical complications.
- Do not use this unit if you are injured.
- Do not use this appliance in the genital area.
- Do not use this unit on your face or head.
- Do not use this unit directly on bone.
- Use the device only on soft and muscle tissue.
- Bruises and contusions can occur even at low frequency and low pressure.
- Check the treated areas after each use and stop immediately if you notice the first signs of pain or discomfort.
- Move the device slowly and not longer than 60 seconds over the respective muscle group.
- Keep fingers, hair or other body parts away from the massage head opening to avoid crushing!
- Do not use the appliance for more than 20 minutes without interruption.
- Let the unit rest for 30 minutes before using it again.
- Keep the appliance away from liquids or heat sources.
- Creams, lotions and oils must only come into contact with the metal attachments provided.
- Keep the ventilation openings free of dirt and dust.
- Do not remove any screws and do not open the unit.
- Only use the original charger supplied.
- Disconnect the charger from the battery after charging and before use.
- Do not use the unit while charging.

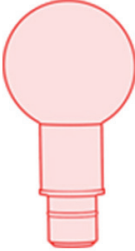
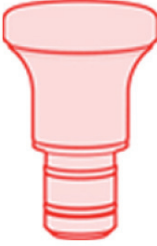
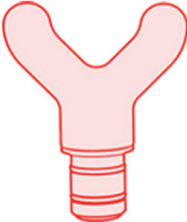
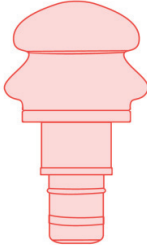
## TECHNICAL DATA

Product code	10038969
Battery capacity	2600 mAh
Power	22-220 W
Frequency	22-55 Hz
Battery life	2-7 hours
Torque	8.5 kg/cm
Amplitude	12 mm

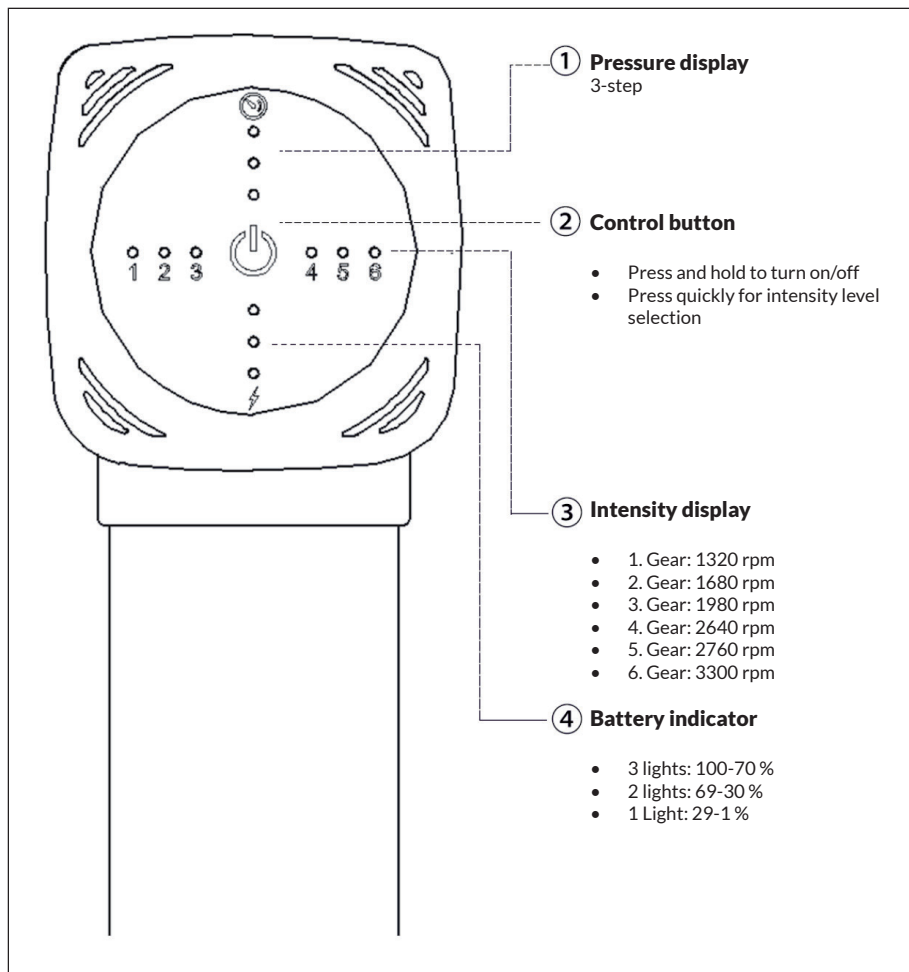
## SCOPE OF DELIVERY

**ESSAYS**

Ball attachment	Flat top	Cone attachment
		
<p>For large muscle groups, such as thighs or glutes.</p>	<p>For general use for all muscle groups.</p>	<p>For the precise treatment of muscles and trigger points.</p>
Fork attachment	Buffer attachment	
		
<p>For medium-sized muscle groups such as the neck, along the spine or along the calves.</p>	<p>For general use for sensitive tissue.</p>	

## CONTROL PANEL AND DISPLAYS

---

## OPERATION

---

- 1 Charge the unit for 5 hours before first use. Do not switch the unit on while charging.
- 2 Select the appropriate massage head and attach it to the massage gun with a twisting motion.
- 3 To remove the attachment, turn it and pull it off.
- 4 Press and hold the control button to switch on the unit.
- 5 Briefly press the control knob to change the intensity level.
- 6 Finally, press and hold the control button again to switch off the device.

---

## DISPOSAL CONSIDERATIONS

---

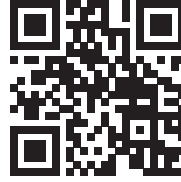
If there is a legal regulation in your country regarding the disposal of electrical and electronic equipment, this symbol on the product or on the packaging indicates that this product must not be disposed of with household waste. Instead, it must be taken to a collection point for the recycling of electrical and electronic equipment. By disposing of this product in accordance with the regulations, you protect the environment and the health of those around you from negative consequences. For information on recycling and disposal of this product, contact your local government or household waste disposal service.

This product contains batteries. If there is a legal regulation in your country regarding the disposal of batteries, the batteries must not be disposed of in household waste. Consult your local regulations for the disposal of batteries. By disposing of this product in accordance with the regulations, you protect the environment and the health of those around you from negative consequences.

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**IMPORTADOR:**

Chal-Tec GmbH, Wallstraße 16, 10179 Berlín, Alemania.

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Chal-Tec UK limited.  
Unit 6 Riverside Business Centre  
Brighton Road  
Shoreham-by-Sea  
BN43 6RE  
Reino Unido

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