CAPITAL SPORTS

ReCrest

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Leg and knee massager

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
Berlin Brands Group UK Ltd
PO Box 1145
Oxford, OX1 9UW
United Kingdom

Dear customer,

Congratulations on the purchase of your device. Please read the following instructions carefully and follow them to prevent potential damage. We accept no liability for damage caused by disregarding the instructions and improper use. Please scan the QR code to access the latest operating instructions and further information about the product.

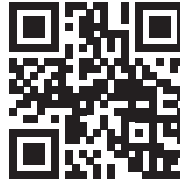CONTENTS

Safety Instructions	12
Device Overview	13
Operation	14
Cleaning and Care	15
Troubleshooting	15
Disposal Considerations	16
Manufacturer & Importer (UK)	16

TECHNICAL DATA

Article number	10040708
Power supply (power adapter)	220-240 V ~ 50/60 Hz
Power consumption	24 W
Preset time	20 minutes


SAFETY INSTRUCTIONS

- Only connect the power supply to sockets that correspond to the voltage specified in the technical data.
- Do not attempt to repair or disassemble the product yourself as there are no user-serviceable parts.
- Do not use the unit in a place where it could get wet or damp.
- Stop using the unit if you feel unwell or if any abnormality occurs.
- Do not use this product in an environment with high temperatures.
- Do not use this product if you are pregnant or menstruating.
- Do not store the unit where it could come into contact with flammable objects, explosive objects or chemicals.
- Do not touch the unit with sharp objects to avoid damaging it.

Do not use this product if any of the following apply to you:

- Your skin is infected, swollen, burnt or otherwise injured.
- You have a pacemaker or other medical implant.
- You suffer from heart disease.
- Your bones are brittle or injured.
- You suffer from cancer.
- You are being treated for your blood pressure or your blood pressure is not normal.
- You suffer from osteoporosis.

DEVICE OVERVIEW

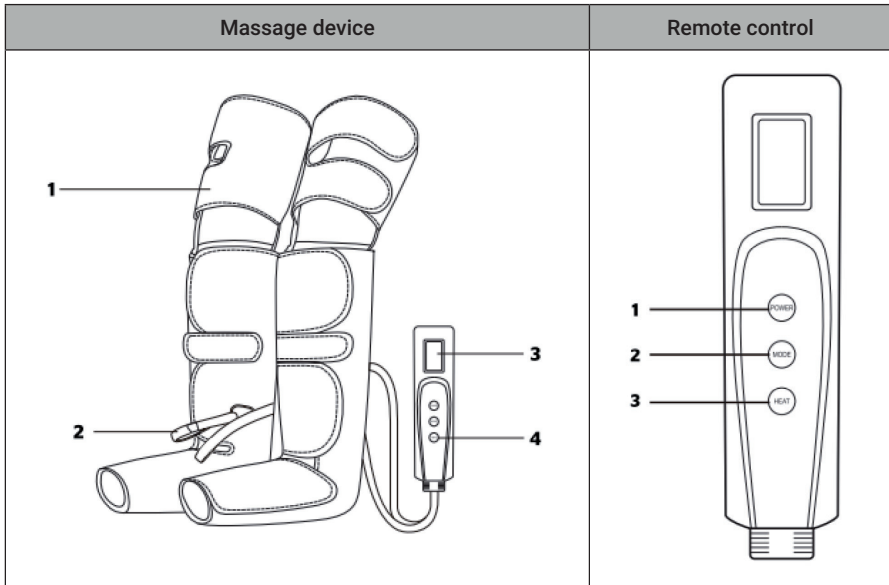Massage device

1	Leg massager	3	LCD display
2	Compressed air hose	4	Control panel

Remote control

1	Switch the device on and off: Press the button briefly to switch the device on. To change the massage level, briefly press the button again. To turn the device off, press and hold the button for 3 seconds.
2	Select mode: Press the button briefly to select one of the following 3 modes: (1) thigh, calf and foot massage, (2) thigh and calf massage only, (3) foot massage only.
3	Select heating function: Press the button briefly to select one of the 3 heating levels.

OPERATION

1. Slip your legs into the leg cuffs and let your legs lie flat on the floor. Note that the leg cuffs are marked L (left leg) and R (right leg).
2. Use the foot strap to secure your foot firmly in the massager.
3. Use the leg strap to fix your leg firmly in the massager.
4. Move your foot to make sure the massager is securely attached.
5. Plug the two compressed air connections into the control unit.
6. Attach the power supply unit to the leg massager and insert the plug into the mains socket.
7. Press the POWER button on the remote control to turn on the device.
8. Then press the POWER button briefly several times to set one of the following intensity levels: LOW-MEDIUM-HIGH.
9. Then press the MODE button to set the desired mode.
10. Then press HEAT button to set the desired heating level.

Note: The appliance switches off automatically after 20 minutes. Do not use the leg massager for more than one hour.

CLEANING AND CARE

- Store the product in a dry place, away from moisture and direct sunlight.
- Wipe the unit with a dry towel after use.

TROUBLESHOOTING

Problem	Possible cause	Solution
The device does not work.	The device is not connected to a power source.	Make sure that the unit is fully connected to the mains.
	The POWER button has not been pressed.	Press the POWER button.
The device suddenly stopped working.	The unit was in operation for more than 20 minutes.	Restart the unit to put it back into operation.

DISPOSAL CONSIDERATIONS

If there is a legal regulation in your country regarding the disposal of electrical and electronic equipment, this symbol on the product or on the packaging indicates that this product must not be disposed of with household waste. Instead, it must be taken to a collection point for the recycling of electrical and electronic equipment. By disposing of this product in accordance with the regulations, you protect the environment and the health of those around you from negative consequences. For information on recycling and disposal of this product, contact your local government or household waste disposal service.

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Chal-Tec GmbH, Wallstraße 16, 10179 Berlino, Germania.