

HALO®

Compact Percussion Massager



Model:239838

Welcome

Thank you for choosing HALO. We design products that help you align, recover, and restore. Use this compact massager for targeted relief. Read the safeguards before first use.

Important Safeguards

1. Read all instructions.
2. This device is not intended to diagnose, treat, cure, or prevent any disease.
3. Consult your physician before use if you experience any pain or have any of the following health conditions:
 - (1) You are pregnant, or have diabetes, retinal problems, or neuropathy.
 - (2) You have internal medical devices in your body.
 - (3) You have other medical conditions or health concerns.
4. This device is charged by a 5V/2A charger or a fast charging power adapter.
5. Fully charge the device before its first use.
6. Do not replace the attachment head while it is running.
7. Do not immerse this device in water or other liquid.
8. Do not use this device to massage your head, chest, genitals, or any area with wounds. Avoid massaging the surgical area during the healing process.


Caution

This device has a built-in rechargeable lithium battery.

1. Do not dispose of the device with household trash.
2. Properly handle it by following your local recycling guidelines.

Specifications


Model: 239838
Dimensions: 95 x 48 x 148mm
Net weight: 390 g
Speed settings: 5
Speed: 5 (1700–3200 rpm, includes Smart Mode)
Battery capacity: 1800 mAh x 2pcs
Auto shutoff: up to 15 min
Power supply: rechargeable lithium battery
Charging time: approx. 200 min
Battery life: up to 6 hours at low speed
Rated voltage: 7.4V
Max. voltage: 8.4V
Max. current: 6.2A
Max. power: 51W
Rated input: 5V---2A
Working noise: 45 dB

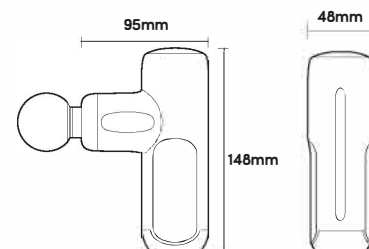

Product Overview

Device (Massager)


- ① Interchangeable attachment head
- ② Grip
- ③ Speed setting indicator
- ④ Battery level indicator
- ⑤ Power on-off / Speed switch button
- ⑥ USB-C port

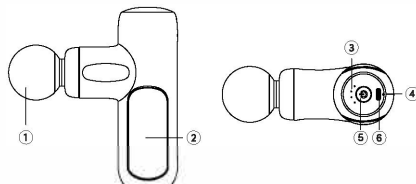

Accessories

1. Attachment heads

- Ball head
- ② Fork head
- ③ Pinpoint head
- ④ Convex head

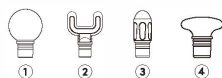

2. USB-C cable



Quick Start Guide

1. Charge the device

Connect the enclosed USB-C cable with a power adapter (not included). The LED display color indicates the battery level.

- (1) The red light flashes - the battery level is low.
- (2) The red light stays on - the battery is in charging.
- (3) The green light stays on - the battery is fully charged.

2. Select an attachment head

- (1) Ball head: Low impact for massaging sensitive and sore areas
- (2) Fork head: Designed for massaging the spine and neck
- (3) Pinpoint head: Pinpoint massage including hands, feet, and knots
- (4) Convex head: Great for massaging large muscle groups and overall use

3. Insert the selected attachment head to the device

Note: Never change an attachment head when the device is on.

4. Turn on the device

Press the Power on-off/Speed switch button ⑤ for 2 seconds to turn on the massager. Initially it will be set at Speed 1.

Note: Press the Power on-off/Speed switch button ⑤ for 2 seconds toggles power on and off.

5. Set your desired speed

After the device is on, tap the Power on-off / Speed switch button ⑤, speed will be shifted up and the next speed light is on. Keep tapping the button ⑤ until your desired speed light is on. The speed setting is done.

Note: The higher the speed, the higher the intensity of the massage. Below is the general guideline for your reference.

- (1) Speed 1 – 1700 rpm: for warming up muscles or very sore muscles
- (2) Speed 2 – 2200 rpm: for soothing tired muscles or sore muscles
- (3) Speed 3 – 2700 rpm: for relaxing stiff muscles
- (4) Speed 4 – 3200 rpm: for massaging very stiff muscles or deep tight knots
- (5) Speed "Smart Mode": four blue lights scrolling display: automatic frequency conversion mode, will automatically adjust the massage strength and speed according to your pressing strength

6. Massage

Apply massage head to your body's soft tissues or target areas.

7. Charging protection

- (1) System protection (low voltage prompt): When the battery is sufficient, a green light will be displayed; When the battery level is low, the red light will flash. Please charge it as soon as possible. If it is not charged, it will automatically shut down for protection later.
- (2) Locked rotor protection: When the locked rotor current reaches the protection value, the motor stops rotating; After entering the locked rotor protection state, short press the function button to resume operation.
- (3) Massage duration protection, the device will automatically shut down after running continuously for 15 minutes (including pause state).

Customer Care

For customer care, contact us at:
customer@halomassagers.com