

Laboratory

Water Distillation

DS 2000

DS 4000

DS 8000

DS 8008

DS 8012

DS 8025

 DragLab

Connect & Go



Laboratory Water Still DS 8008

Professional, fully automatic water distiller with a production capacity of 8 liters per hour and an integrated 16-liter storage tank. Constructed from high-quality stainless steel, this distiller is designed to produce high-purity distillate autonomously, eliminating the need for continuous personnel supervision.

Effortless Water Purity With DragLab Distillation Solutions

DragLab Distillation Units are the most effective and reliable devices for producing pure water, requiring only a connection to raw water and a power source. This process involves transforming raw water (tap water) into steam and then back into liquid, effectively separating water from impurities with higher boiling points. The resulting steam, once condensed, produces distilled water with a purity of approximately 99.5%, removing salts, organic substances, microorganisms, pyrogens, and bacteria. DragLab water distillation units require minimal maintenance needing only regular cleaning. When carbon dioxide from the air dissolves in distilled water, its pH becomes slightly acidic, 5 to 5.6. Distilled water should be used promptly or stored in airtight conditions to maintain purity.



Application:

- Biotechnology, cell culture and media preparation
- Testing, analysis and research laboratories
- Microbiology and Bacteriology sample preparation
- Food and beverage production
- Sterilization, cleaning and Medical equipment.
- Cosmetics, skincare products production.
- Electronics and Automotive.
- Pharmaceutical manufacturing and drug formulation.

A vertical bar on the left side of the page, divided into three colored segments: black at the top, red in the middle, and yellow at the bottom.

Quality and Expertise



DragLab is certified to ISO 9001:2015, reflecting our unwavering commitment to quality and excellence. We continuously optimize our production processes, ensuring that every step is efficient and effective. Through ongoing development and innovation, we strive to meet and exceed industry standards. Our customer-centric approach prioritizes the needs and satisfaction of our clients, driving us to maintain rigorous quality assurance protocols. This dedication ensures that laboratory equipment manufactured by DragLab consistently meets the highest standards of quality and reliability.

DragLab Laboratory Devices are meticulously designed and manufactured from high-quality materials to meet the most rigorous requirements of our customers. Our commitment to excellence and innovation drives continuous development and improvement, informed by daily usage and ongoing technological advancements. 'Quality and Expertise' is not just about creating impressive product designs; it reflects our company's policy of prioritizing high preparedness and customer-centric service as our fundamental objective. This dedication extends to all DragLab products, which comply with current European standards and bear the CE mark.




Laboratory Water Still DS 8008

Professional, fully automatic water distiller with a production capacity of 8 liters per hour and an integrated 16-liter storage tank. Constructed from high-quality stainless steel, this distiller is designed to produce high-purity distillate autonomously, eliminating the need for continuous personnel supervision.



Dimensions:

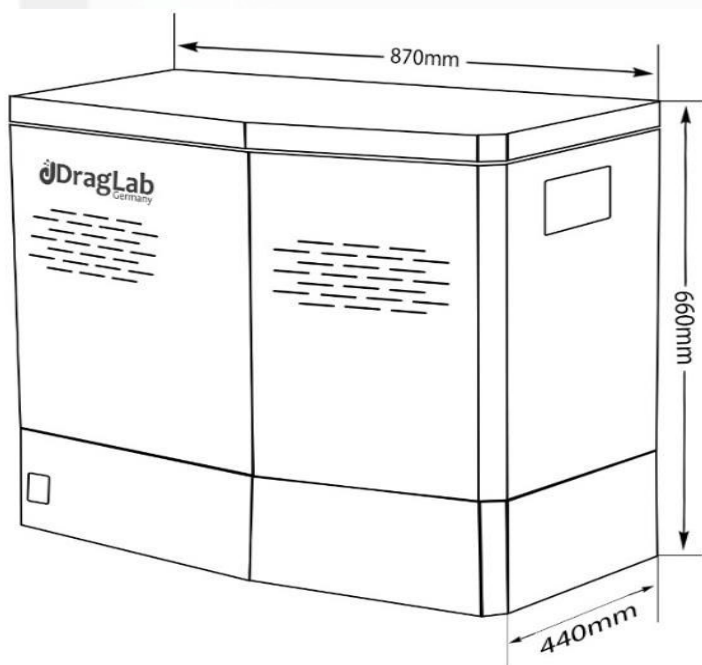

Features:

Distillation Capacity:

Produces 8 L/h of distilled water.

Storage Tank:

Integrated with an 16-liter storage tank.

Purity:

High-purity distilled water with conductivity of approximately 2.5 $\mu\text{S}/\text{cm}$ at 25 °C.

Water Quality:

Pyrogen-free, low-gas, and bacteria-free water.

Protection:

Automatic thermostatic low-water cut-off to protect the heating element in case of low water.

Operation:

Automatic water switching for continuous functionality.

Energy Efficiency:

Energy and water-saving through an automatic distillation process.

Display:

Screen indicating water level in the storage tank, heater operation status, water supply, error status check, and user instructions.

Maintenance:

Easily accessed evaporator tank for effortless cleaning and maintenance.

Material:

All water contact materials are made of stainless steel (material no.: EN 1.4301).

Heating Element:

Made of stainless steel (material no.: EN 1.4301).

Housing:

Galvanized steel with electrostatic powder coating.

Installation:

Suitable for both bench and wall mounting.

Accessories:

All parts and accessories for installation included.

Laboratory Water Still DS 8008

Specifications:

Capacity

Production Capacity	8 liter / hour
Storage Tank	Integrated 16 Litter storage tank
Consumption	80 liter / hour
Conductivity	Approximately 2,5 µS/cm at 25 °C

Control technology

Fully Automatic	Automatic thermostatic power cut-off to protect the heating element in case of low water
Fully Automatic	Automatic power activation when the water level reaches the safety line in the evaporator tank to protect the heating element.
Fully Automatic	Automatic refilling of the storage tank after distillate withdrawal , electronically controlled.
Water supply	Built-in solenoid valve connection for the evaporator and cooler, with automatic switching on/off based on the water level in the evaporator tank.

Standard equipment

Internal material	The tanks (Evaporating tank & Storage tank) are made of stainless steel (material no. EN 1.4301)
Internal material	The condenser (cooler) is made of stainless steel (material no. EN 1.4876).
Internal material	The heating element is made of stainless steel (material no. EN 1.4876).
External material	The housing is made from galvanized steel and coated with electrostatic powder.
Degassing system	The unit includes a top vent for effective degassing, ensuring pyrogen-free, low -gas, and bacteria-free water.
Accessing interior	The evaporator and storage tanks are easily accessible for cleaning and maintenance, facilitating efficient cleaning procedures.
Installation	All necessary parts and accessories for installation are included.
Distillate water withdrawal	Distillate water can be withdrawn by continuously or push mode from front tap.
Main switch	The main power switch button located on the front of the unit
Screen	The device features a screen indicating water level in the storage tank, heater operation status, water supply, error status check, and user instructions.
Mounting	Suitable for both benchtop and wall mounting.

Safety

Water	Low water level sensor with automatic power cut-off to protect the heating element.
Heating	Overheating cut-off
Screen	The screen displays startup tests for the heating element and provides visual alarms for faults and errors.
Degassing	pyrogen-free, low-gas, and bacteria-free water.
Power connection	Power cable with 5 pins three-phase CEE plug.

Laboratory Water Still DS 8008

Specifications:

Connection

Water inlet	The cooling water inlet is located on the right-hand side of the unit, with a diameter of 3/4 inch (Ø 19.05 mm).
Water pressure range	From 2 bars to 10 bars
Water outlet	The cooling water outlet is located on the right-hand side of the unit, with a diameter of 2/5 inch (inner Ø 10 mm).
Drainage tap	The drainage tap is located on the right-hand side of the unit, with a diameter of 2/5 inch (inner Ø 10 mm).
Distillate water	Distillate water can be withdrawn by continuously or push mode from front tap.
Power	The main connection cable on the left-hand side of the unit

Electrical data

Voltage	380 volts, 50/60Hz
Watt	6000 watts

Shipping information

Customs tariff number	8419 4090
Packaging Dimension	94 W X 54 D X 84 H (CM)
Net weight	40 Kg
Gross weight carton	46 Kg

Order Information:

Description	Model	Article number
Water Still – Capacity: 8 L/H, 16 L storage Tank – 380V, 50/60Hz	DS 8008	5008.000

