


Barcode Plus™ Slim|Handwritten

The Barcode Plus™ Slim Handwritten Blood Band has plenty of space for handwriting patient information and includes a shield to protect patient information. The large Blood Bank ID (BBID) makes it quick and easy to read, and the barcodes allow greater efficiency when entering the BBID in your Laboratory Information System (LIS).

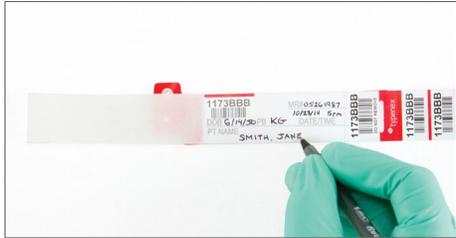

1. **VERIFY** the patient's identity.

WRITE the patient's information over the prompts on the band using a ballpoint pen.

NOTE: Do not apply shield over handwritten information at this time.



2. **WRAP** the band tail around the patient's extremity to size.


SNAP the button closure to secure.

CUT the excess band tail material with scissors.



3. **DRAW** the pre-transfusion specimen per your protocol.

PEEL the handwritten patient information label from the band.

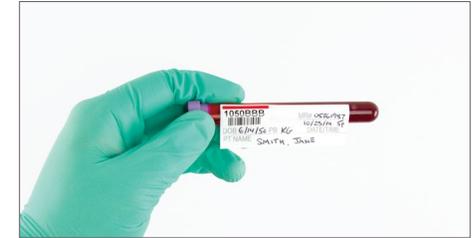

4. **PLACE** the handwritten patient information label on the specimen tube lengthwise.



5. **PEEL** the liner from the shield.

COVER the copy of the handwritten patient information with the shield.

SEAL the edges down using your fingers.



6. **LABEL** any forms at the patient's bedside as required by your protocol.

REMOVE the tail of stickers from the band at the perforation.

ATTACH the tail of stickers to the specimen tube, if desired.



7. **PEEL** the adhesive tab on the underside of the extra stickers.

SECURE the extra stickers to the band.



8. **TRANSPORT** the specimen to the blood bank.

SCAN the BBID into the LIS or enter it manually.

Visit www.typenexmedical.com or call us at 866.897.3639 to learn more.