

ProteoSpin™ Inclusion Body Protein Isolation Maxi Kit

The ProteoSpin™ Inclusion Body Isolation Maxi Kit facilitates the isolation of recombinant proteins in the form of inclusion bodies from *E. coli*. The kit includes reagents specially formulated to achieve rapid and high-quality purification of inclusion body proteins using three processes:

- Lysis of bacterial cells to release inclusion bodies in solid form
- Solubilization of purified inclusion bodies
- Purification of the recombinant protein using spin column chromatography



With optimized reagents and streamlined processes, the ProteoSpin™ Inclusion Body Isolation Maxi Kit significantly reduces time and labour for isolating and purifying proteins of interest in a larger-scale format. The ProteoSpin™ Inclusion Body Isolation Maxi Kit employs spin-column chromatography using Norgen's proprietary resin as an ion-exchanger. Each spin column is able to purify up to 12 mg of recombinant proteins from 100 mL of culture. The kit is designed to purify both acidic and basic proteins.

Kit Specifications			
Maximum Culture Volume	100 mL	Minimum Elution Volume	4 mL
Yield from 100 mL of Culture	Up to 12 mg	Time to Process 1 Sample	1 hour

ProteoSpin™ Inclusion Body Protein Isolation Maxi Kit Benefits

Three kits in one for acidic and basic proteins	The kit provides essential reagents for cell disruption, inclusion body solubilization, and purification using spin column chromatography. The kit includes solutions and protocols for use with either acidic or basic proteins.
Gentle disruption procedure provides high purity final product	Cell lysis is accomplished through non-ionic detergent chemical disruption. The use of a needle and syringe reduces viscosity and facilitates purification, producing proteins that are often >95% pure. Proteins are then ready for SDS-PAGE, 2D gels, mass spectrometry analysis or other procedures.
Convenient process volume	Analyze 100 mL of microbial cell culture yielding up to 12 mg of protein.
Fast and easy processing	Process time is about 1 hour for a single sample.
Long kit shelf life	The kit uses an enzyme-free process with stable solutions for longer shelf life.
Suitable for downstream applications	Final elution performed with volatile or non-volatile buffers suitable for: <ul style="list-style-type: none"> • Mass spectrometry • SDS-PAGE • Refolding experiments

ProteoSpin™ Inclusion Body Protein Isolation Maxi Kit

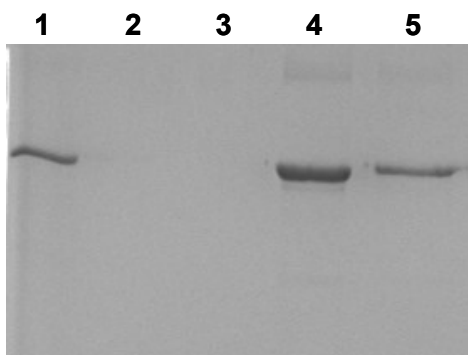


Figure 1. No Loss of Proteins when Isolating a Basic Protein

No proteins are lost during the isolation of a recombinant 30 kDa basic protein using the ProteoSpin™ Inclusion Body Protein Isolation Maxi Kit. Lane 1 is the input, Lane 2 is the flowthrough from the input, Lane 3 is the wash, Lane 4 is the first elution and Lane 5 is the second elution.

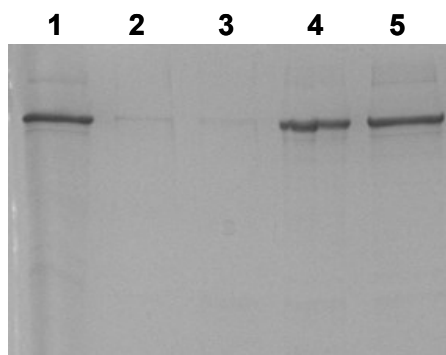


Figure 2. No Loss of Proteins when Isolating An Acidic Protein

No proteins are lost during the isolation of a recombinant 70 kDa acidic protein using the ProteoSpin™ Inclusion Body Protein Isolation Maxi Kit. Lane 1 is the input, Lane 2 is the flowthrough from the input, Lane 3 is the wash, Lane 4 is the first elution and Lane 5 is the second elution.

Inclusion Body Protein Isolation Maxi Kit Contents

1. Wash Solution C
2. Wash Solution N
3. Binding Buffer A
4. Binding Buffer N
5. Elution Buffer C
6. Protein Neutralizer
7. Cell Lysis Reagent
8. IB Solubilization Reagent
9. Syringes (10 mL, Luer-Lok™ Tip)
10. Needles (18G x 1½ inch)
11. Maxi Spin Columns (filled with SiC) inserted into 50 mL collection tubes
12. Final Elution Tubes, 50 mL
13. Product Insert

Shipping Conditions

The ProteoSpin™ Inclusion Body Protein Isolation Maxi Kit is shipped at room temperature.

Storage Conditions

The Cell Lysis Reagent and the IB Solubilization Reagent should be stored at 4°C upon receipt of the kit. All other unopened solutions should be stored at room temperature. Once opened, the solutions should be stored at 4°C when not in use, except for the Basic and Acidic Binding Buffers, which should be stored at room temperature.

Customer-Supplied Reagents and Equipment

- Benchtop centrifuge (capable of spinning 50 mL conical tubes)
- Centrifuge (capable of spinning 40 mL tubes and 250 mL bottles)
- Centrifuge tubes (40 mL)
- 50 mL conical tubes
- Sterile 100 mL bottles
- Sterile, deionized water or Milli-Q® water

Cat #	Description	Quantity
17700	ProteoSpin™ Inclusion Body Isolation Maxi Kit	4 preps