
502 Elastic Ankle support

DE	Gebrauchsanweisung	4
EN	Instructions for use	8
FR	Instructions d'utilisation	12
IT	Istruzioni per l'uso	16
ES	Instrucciones de uso	20
NL	Gebruiksaanwijzing	24
SV	Bruksanvisning	28
DA	Brugsanvisning	31
NO	Bruksanvisning	35
FI	Käyttöohje	39
PL	Instrukcja użytkowania	43
ET	Kasutusjuhend	47
RU	Руководство по применению	51

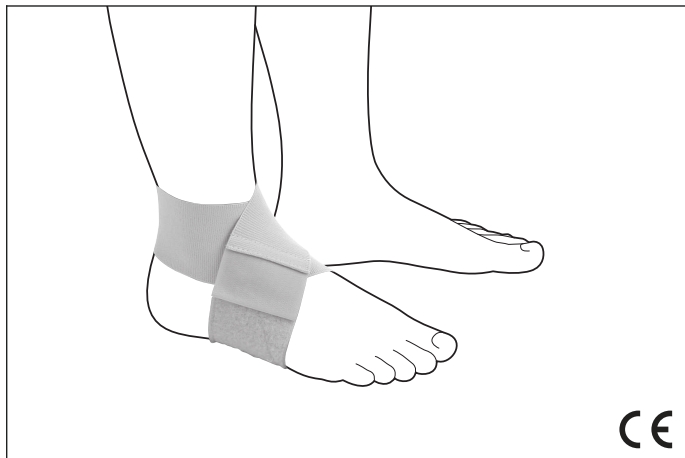	Artikelnummer / Article number	Größe / Size	Umfang / Circumference	
			cm	inch
	502=XS	XS	49-51	19.3-20.1
	502=S	S	51-53	20.1-20.9
	502=M	M	53-55	20.9-21.7
	502=L	L	55-57	21.7-22.4
	502=XL	XL	57-59	22.4-23.2
	502=XXL	XXL	59-61	23.2-24.0
Seite / Side: beidseitig / suitable for either side				

Material	PES, PA, EA
-----------------	-------------

INFORMATION

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INFORMATION

Date of last update: 2026-09-02

- ▶ Please read this document carefully before using the product and observe the safety notices.
- ▶ Instruct the user in the safe use of the product.
- ▶ Please contact the manufacturer if you have questions about the product or in case of problems.
- ▶ Report each serious incident related to the product to the manufacturer and to the relevant authority in your country. This is particularly important when there is a decline in the health state.
- ▶ Please keep this document for your records.

These instructions for use provide important information on the application of the 502 elastic ankle support.

2 Intended use

2.1 Indications for use

The support is intended **exclusively** for fittings of the lower limbs and **exclusively** for contact with intact skin.

The support must be used according to the indication.

2.2 Indications

- Feeling of instability in the ankle, e.g. in case of rheumatic complaints or following a sprain.
- Irritation and swelling of the ankle joint

Indications must be determined by the physician.

2.3 Contraindications

2.3.1 Absolute Contraindications

Not known.

2.3.2 Relative Contraindications

The following indications require consultation with a physician: skin diseases/injuries, inflammation, prominent scars that are swollen, reddening

and hyperthermia of the fitted limb; pronounced varicose veins, especially with return flow impairment, lymphatic flow disorders, including unclear soft tissue swelling distant to the body area to which the support will be applied; sensory and circulatory disorders in the foot area, e.g. in case of diabetic neuropathy.

2.4 Lifetime

The product is designed for a maximum lifetime of **6 months**.

2.5 Effects

The support stabilises the ankle and can relieve pain.

3 Safety

3.1 Explanation of warning symbols

 CAUTION	Warning regarding possible risks of accident or injury.
---	---

NOTICE	Warning regarding possible technical damage.
---------------	--

3.2 General safety instructions

CAUTION

Contact with heat, embers or fire

Risk of injury (such as burns) and risk of product damage.

- ▶ Keep the product away from open flames, embers and other sources of heat.

CAUTION

Reuse on other persons and improper cleaning

Skin irritation, formation of eczema or infections due to contamination with germs

- ▶ The product may be used by one person only.
- ▶ Clean the product regularly.

NOTICE

Improper use and changes

Change in or loss of functionality as well as damage to the product

- ▶ Use this product with care and only for its intended purpose.
- ▶ Do not make any improper changes to the product.

NOTICE

Contact with oils, salves, lotions or other products that contain oils or acids

Insufficient stabilisation due to loss of material functionality

- ▶ Do not expose the product to oils, salves, lotions or other products that contain oils or acids.

4 Handling

INFORMATION

- ▶ The daily duration of use and period of application are generally determined by the physician.
- ▶ The initial fitting and application of the product must be carried out by qualified personnel.
- ▶ Instruct the patient in the handling and care of the product.
- ▶ Instruct the patient to see a physician immediately if any exceptional changes are noted (e.g. worsening of the complaint).

4.1 Size Selection

- 1) Use a tape measure to make an "8" around the ankle and midfoot (size chart).
- 2) Measure the ankle circumference above the lateral malleolus and the midfoot circumference (metatarsus) to determine the overall circumference.
- 3) Select the size of the support (see size chart).

4.2 Application

⚠ CAUTION

Incorrect or excessively tight application

Skin irritation due to overheating, local pressure points due to tight application to the body

- ▶ Ensure that the product is applied properly and fits correctly.
- ▶ When in doubt, do not continue to use the product if skin irritation occurs.

- > Open all hook-and-loop closures on the support.
 - > The patient is sitting on the front edge of a chair.
- 1) Apply the middle of the support to the Achilles tendon and cross the two ends over the arch of the foot in front (see fig. 1).

- 2) Pull the end of the support with the soft hook-and-loop closure under the foot to the outer side.
- 3) Close the end of the support with the rough hook-and-loop closure together with the soft hook-and-loop closure on the outside of the arch (see fig. 2).

4.3 Cleaning

NOTICE

Use of improper cleaning agents

Damage to the product due to use of improper cleaning agents

- Only clean the product with the approved cleaning agents.

Clean the support regularly:

- 1) Fasten all hook-and-loop closures.
- 2) Recommendation: Use a laundry bag or net.
- 3) Wash the support in warm water at **40 °C** with standard mild detergent. Do not use fabric softener. Rinse thoroughly.
- 4) Allow to air dry. Do not expose to direct heat sources (e.g. sunlight, stove or radiator).

5 Disposal

Dispose of the product in accordance with national regulations.

6 Legal information

All legal conditions are subject to the respective national laws of the country of use and may vary accordingly.

6.1 Liability

The manufacturer will only assume liability if the product is used in accordance with the descriptions and instructions provided in this document. The manufacturer will not assume liability for damage caused by disregarding the information in this document, particularly due to improper use or unauthorised modification of the product.

6.2 CE conformity

The product meets the requirements of Regulation (EU) 2017/745 on medical devices. The CE declaration of conformity can be downloaded from the manufacturer's website.

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