

GL 48 mg/dL
& mmol/L

beurer
medical



HealthManager



test strip
ejection



target
scale



Coloured target scale and display to help you to interpret the measurement



Electronic diary available on PC

Backlit display with self-explanatory user menu

Easy to use: handy test strips, easy measurement marking for optimised diabetes management: before eating, after eating, general marking

Blood volume check: measurement only possible with sufficient sample of blood

Starter set: GL 48 measuring device, LD 04 lancing device, 10 lancet needles, 10 test strips, USB cable, case, instructions for use, brief instructions, blood glucose diary

Ketone warning on display

GlucoMemory software in German and English: simple overview of results and trends in target range

„beurer HealthManager“ PC software for download free of charge at: www.beurer.com > SERVICE > Download

Suitable for self-testing

Required volume of blood: 0.5 microlitres

Duration of measurement: approx. 5 seconds

Code-free

Measured value memory for 900 measurements with date, time and measurement marking

Calculating the average value for 7, 14, 30 or 90 days for
- all measurements - marked measurements

Blood sample: fresh whole blood (capillary) or venous whole blood

Plasma-calibrated

Measurement range: 20 – 630 mg/dL / 1,1 – 35,0 mmol/L

Fully fulfils EN ISO 15197:2015 (compulsory since May 2016)

Automatic switch-off

Product dimensions (L/H/D): 95 x 53 x 15 mm

Product weight: 52 g

Button cell batteries (CR2032)




GL 48 mg/dL
& mmol/L

beurer
medical



CONNECTION

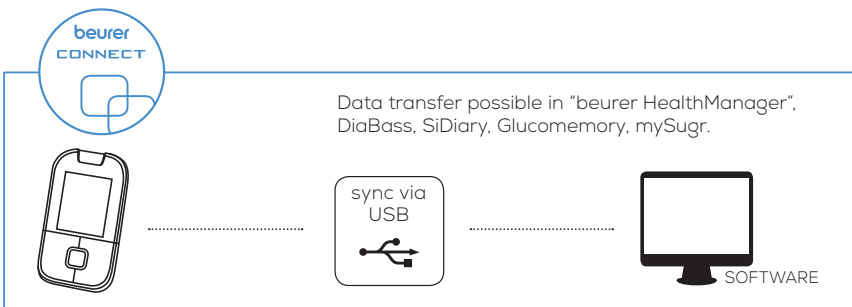

GL48_2020-06-15_03_DS_EN_BEU Subject to errors and changes

* tested online application "beurer HealthManager" Cloud