


Vascular Intervention // **Coronary**  
Covered Coronary Stent System

# PK Papyrus



Exceptional deliverability<sup>1,2</sup>

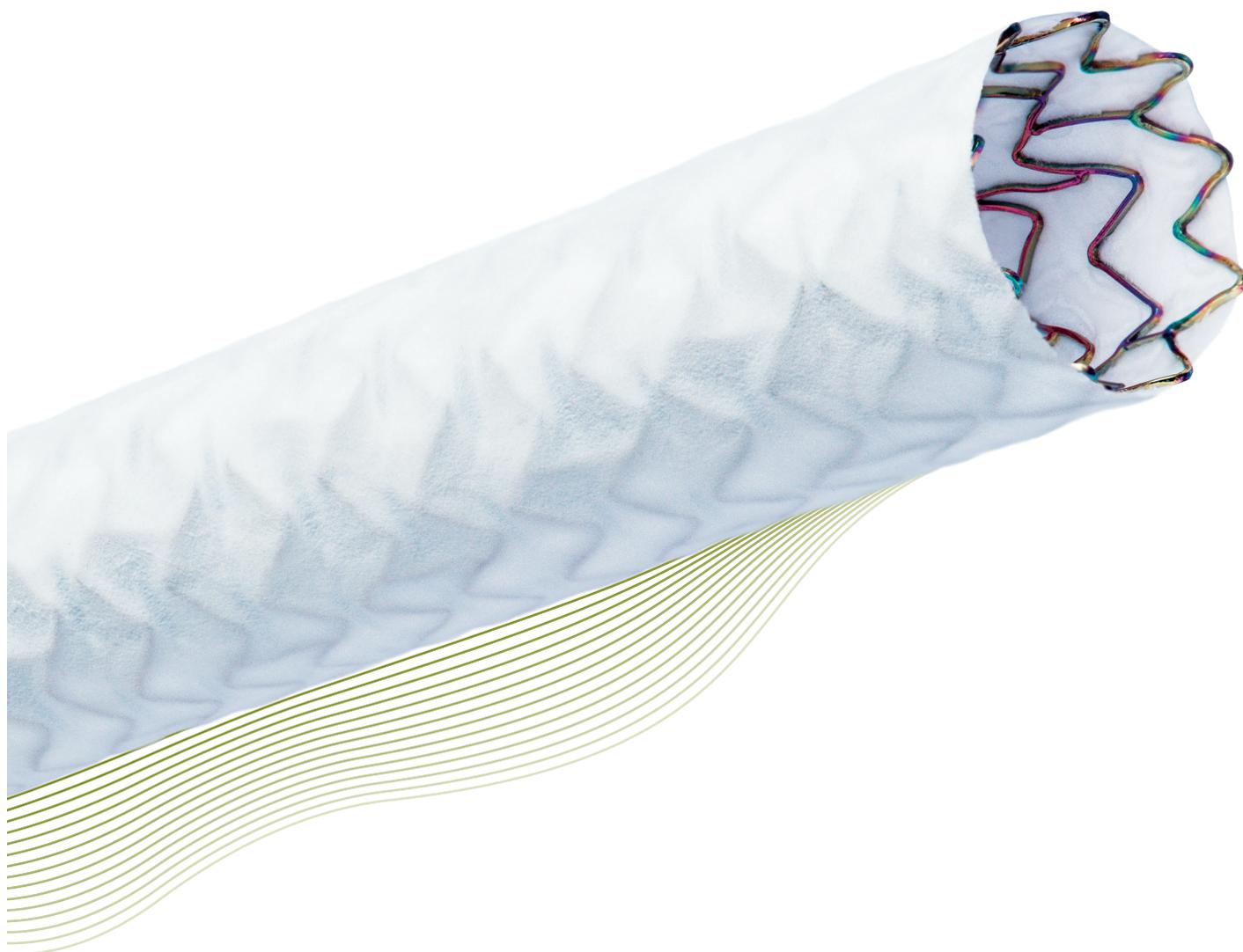

Covered single stent design



Designed to save lives when  
seconds count<sup>3</sup>

Humanitarian Device. Authorized by Federal law for use in the treatment of acute perforations of native coronary arteries and coronary bypass grafts in vessels 2.5 to 5.0 mm in diameter. The effectiveness of this device for this use has not been demonstrated.

 **BIOTRONIK**  
excellence for life

# PK Papyrus

Designed to deliver more like a conventional stent<sup>1, 2</sup>

## Superior design for exceptional deliverability<sup>1, 2</sup>

### Lowest crossing profile<sup>1</sup>


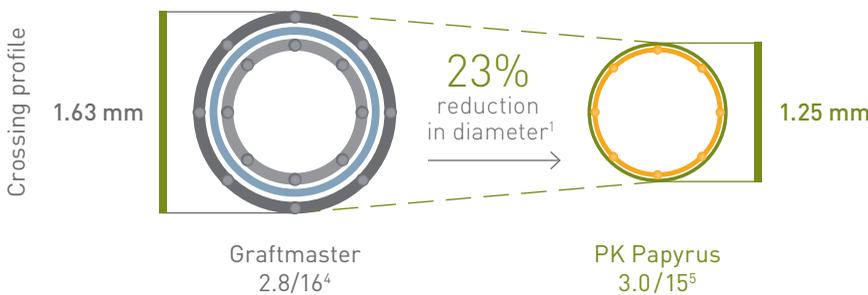



Illustration depicts crimped devices prior to inflation

### Superior flexibility<sup>2</sup>

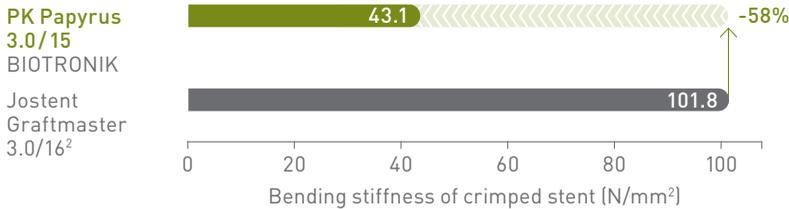


## Covered single stent design


- BIOTRONIK's ultrathin strut stent platform [Cobalt Chromium].
- Highly elastic membrane capable of sealing coronary artery perforations.<sup>6</sup>

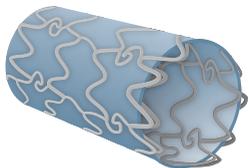


Jostent Graftmaster  
Traditional sandwich stent design

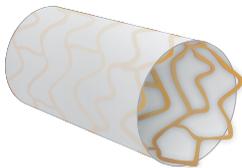

PK Papyrus  
Covered single stent design

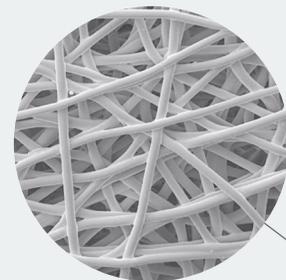

Polyurethane fibers  
(2000x) → 10 μm

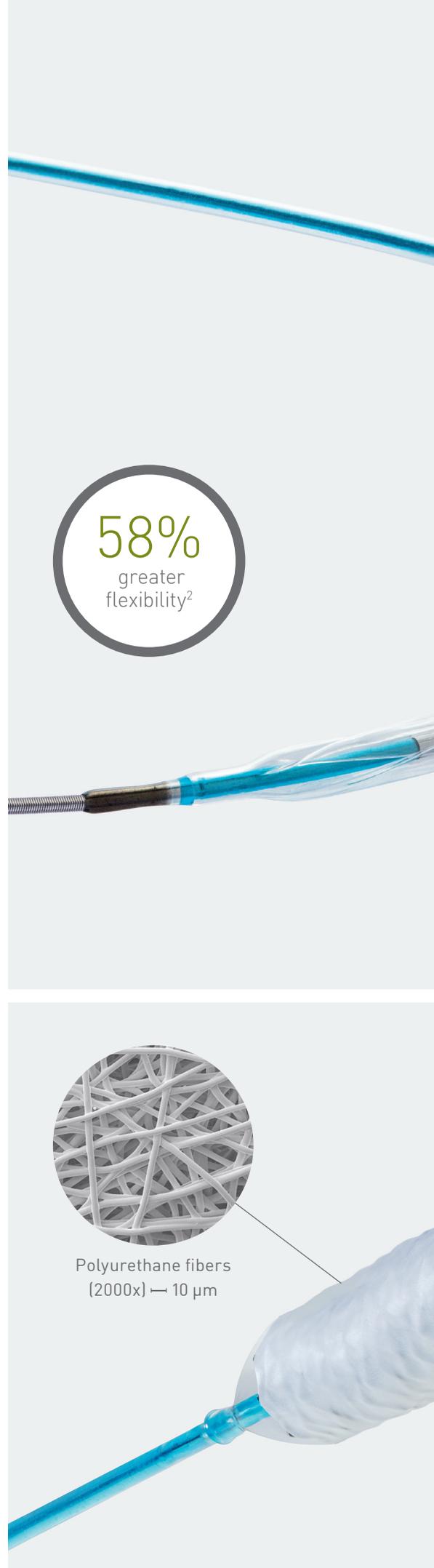

The only  
**5F**  
compatibility\*



**47%**  
less time to  
deliver<sup>6</sup>



# Time to expand your options with **PK Papyrus**

## Broadest range of sizes on the US market<sup>7</sup>

- The only 5F compatible\* covered coronary stent with the broadest range of sizes.<sup>7</sup> For main sizes – no need for guide catheter upgrade.\*
- First FDA approved 2.5 mm diameter.

## Designed to save lives when seconds count<sup>3</sup>

### Shorter time to deliver<sup>6</sup>

Single center, retrospective investigation of 61 patients treated with covered coronary stents.<sup>8</sup>



PK Papyrus  
n = 22



Jostent Graftmaster  
n = 39

Median time to deliver (p=0.001)

\*5F compatible for  $\varnothing$  2.5-4.0 mm; 6F compatible for  $\varnothing$  4.5-5.0 mm

# PK Papyrus

Vascular  
Intervention  
Coronary



Indicated for the treatment of acute perforations of native coronary arteries and coronary bypass grafts in vessels 2.5 to 5.0 mm in diameter\*

Technical Data	Stent
Stent cover material	Non-woven, electrospun polyurethane
Stent cover thickness	90 µm
Stent material	Cobalt chromium (L-605) with <b>proBIO</b> (amorphous silicon carbide) coating
Maximum stent expansion diameter	ø 2.5 - 3.0 mm: 3.50 mm; ø 3.5 - 4.0 mm: 4.65 mm; ø 4.5 - 5.0 mm: 5.63 mm
Delivery system	
Guide wire diameter	0.014"
Usable catheter length	140 cm
Recommended guide catheter	ø 2.5 - 4.0 mm: 5F (min. I.D.** 0.056"); ø 4.5 - 5.0 mm: 6F (min. I.D.** 0.070")
Nominal pressure (NP)	ø 2.5 - 3.5 mm: 8 atm; ø 4.0 - 5.0 mm: 7 atm
Rated burst pressure (RBP)	ø 2.5 - 4.0 mm: 16 atm; ø 4.5 - 5.0 mm: 14 atm

\*\*I.D. = Inner Diameter

Compliance Chart	Inflation pressure	Stent inner diameter (mm)					
		2.5	3.0	3.5	4.0	4.5	5.0
Nominal pressure (NP)	7	-	-	-	4.01	4.55	4.93
Nominal pressure (NP)	8	2.52	2.99	3.53	4.14	4.69	5.09
	9	2.59	3.07	3.63	4.26	4.82	5.23
	10	2.65	3.15	3.71	4.35	4.91	5.34
	11	2.70	3.21	3.77	4.43	4.99	5.43
	12	2.74	3.26	3.82	4.49	5.06	5.50
	13	2.77	3.30	3.86	4.54	5.11	5.56
Rated burst pressure (RBP)	14	2.80	3.34	3.90	4.59	5.16	5.61
	15	2.83	3.37	3.93	4.63	-	-
Rated burst pressure (RBP)	16	2.86	3.40	3.96	4.67	-	-

Ordering Information	Stent ø (mm)	Catheter length 140 cm Stent length (mm)		
		15	20	26
5F	2.5	434887	434893	-
	3.0	434888	434894	434899
	3.5	434889	434895	434900
	4.0	434890	434896	434901
6F	4.5	434891	434897	434902
	5.0	434892	434898	434903

1. Compared to Graftmaster 2.8/16 (BIOTRONIK data on file); 2. Compared to Jostent Graftmaster 3.0/16 (BIOTRONIK data on file); 3. Broad range of sizes available on the US market; 4. Data obtained from Graftmaster Coronary Stent Graft System Brochure 11/13/12; 5. PK Papyrus 3.0/15 (BIOTRONIK data on file); 6. Hernandez-Enriquez M, et al. Outcomes after use of covered stents to treat coronary artery perforations. Comparison of old and new-generation covered stents. J Interv Cardiol. 2018; 1-7; 7. Compared to Graftmaster based on the broader range of sizes available on the US market; 8. Population is representative of real world interventional practice and was not a randomized prospective clinical trial.

Jostent and Graftmaster are registered trademarks of the Abbott Group of Companies.

\*Indication as per IFU.

Information on devices manufactured at companies other than BIOTRONIK was gathered from multiple sources. However, it has not been verified by the vendors and we cannot guarantee its accuracy.

**Manufacturer:**  
BIOTRONIK AG  
Ackerstrasse 6  
8180 Bülach, Switzerland  
Tel +41 (0) 44 8645111  
Fax +41 (0) 44 8645005  
info.vi@biotronik.com  
www.biotronik.com

**Distributor:**  
BIOTRONIK, Inc.  
6024 Jean Road  
Lake Oswego, OR 97035-5369  
Tel (800) 547-0394 (24-hour)  
Fax (800) 291-0470  
www.biotronik.com

442278/A/US/Nov\_2018

© 2018 BIOTRONIK AG – All rights reserved.  
Specifications are subject to modification, revision and improvement.

**BIOTRONIK**  
excellence for life