

ROOT CANAL SHAPING



MicroMega

2Shape mini


Rational preparation approach

micromega.coltene.com

 **COLTENE**

- Preserving pericervical dentin
- Removing minimum of tooth structure¹
- Respecting the root canal anatomy²
- Creating optimal space for a deep disinfection and a proper obturation


Protocol*

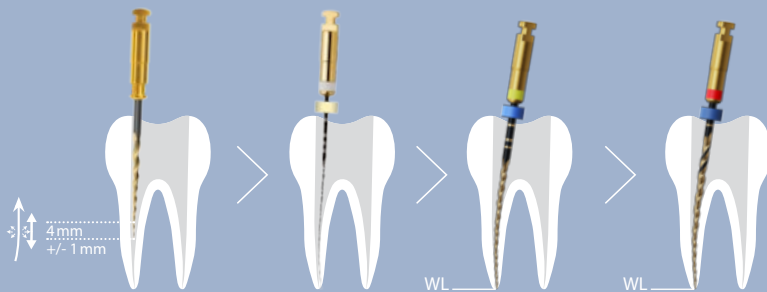

MicroMega:	One Flare	One G	2Shape TS1 mini	2Shape TS2 mini
Speed of rotation:	250-400 rpm	250-400 rpm	250-400 rpm	250-400 rpm
Max. torque:	3 N.cm	1.2 N.cm	2 N.cm	2.5 N.cm

Complete your treatment procedure with matching Guttapercha and Paper Points.

OPTIONAL FINISHING

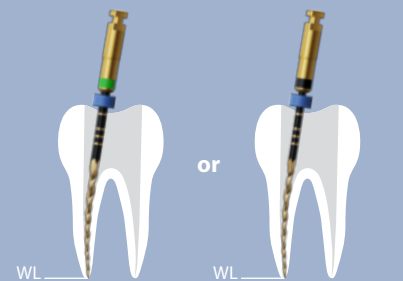

2Shape F35 mini	2Shape F40 mini
250-400 rpm	250-400 rpm
2.5 N.cm	2 N.cm

Features

- Small size/taper
- Ø 1 mm wire
- T.Wire heat treatment
- Cutting Efficiency
- 2 files to shape
- Sterile blister
- 3 lengths: 21, 25, 31 mm
- Multiple use



For more information about MicroMega 2Shape mini, please scan the QR code.

Order Information

MicroMega 2Shape mini, 5 pcs/pack (B5)

20953211	B5 2Shape TS1 mini N25 4% L25 Classics S
20953212	B5 2Shape TS2 mini N25 6% L25 Classics S
20953213	B5 2Shape F35 mini N35 6% L25 Classics S
20953214	B5 2Shape F40 mini N40 4% L25 Classics S

(1) Survey analysis "Minimally Invasive", Caroline Dort, 2022-01-13

(2) PCD2, Shaping Ability of Two-shape and ProTaper Gold Files by using Cone-beam Computed Tomography, The Journal of Contemporary Dental Practice, March 2019;20(3):330-334. 1 Shalini Singh, 2 Tarun Gupta, 3 Vinisha Pandey, 4 Harshit Singhania, 5 Pooja Pandey, 6 Shweta Gangavane

* Extract from the instructions for use.

© COLTENE – micromega.coltene.com

Coltene/Whaledent Ltd.
The President Suite
Kendal House, Victoria Way
Burgess Hill, West Sussex
RH15 9NF/GB
T +44 1444 235486
info.uk@coltene.com

Coltene/Whaledent AG
Feldwiesenstrasse 20
CH-9450 Altstätten
T +41 71 757 5300
info.ch@coltene.com

Micro-Mega SA
12 rue du Tunnel
FR-25000 Besançon
T +33 3 81 54 42 42
commercial@micro-mega.com

COLTENE