



ZARC

ZEvo mini

Endodontic motor with continuous
and reciprocating rotation

**ENDODONTICS MADE EVEN EASIER,
AND IN A MINI SIZE**



FROM ODONTOLOGISTS
TO ODONTOLOGISTS

zEvo mini

Z-Evo mini performs root canal treatments using endodontic instruments in continuous rotation controlled by torque and reciprocating motion.




360° ADJUSTABLE HEAD

Its 360° adjustable head allows easy access to any canal.



SMALLER SIZE

Just as powerful as the classic Z-Evo but in a smaller size.

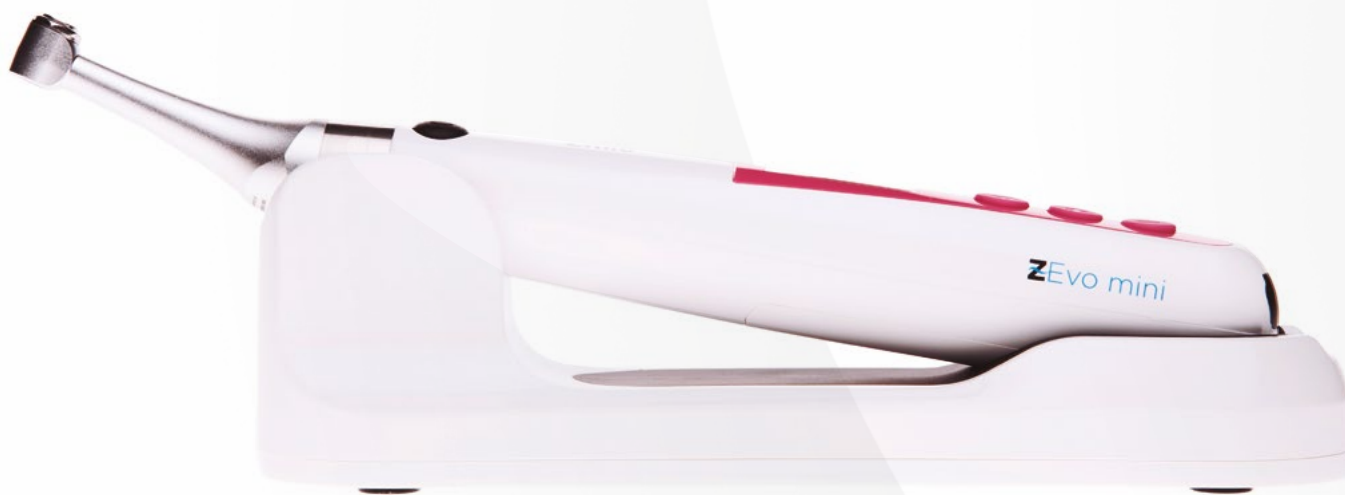

LONG-LASTING BATTERY

It has a long-lasting 800 mAh lithium battery that allows the device to be used while charging. The convenient charging base included makes for an easy charging process.




SAFE TORQUE CONTROL

The immediate torque during use is displayed on the screen, so the dentist can regulate the parameters and adapt them to each system. This avoids exceeding maximum limits and file breakage.



RECIPROCATING MOTOR

Z-Evo mini performs rotary and reciprocating movements, it adapts to any file system on the market.

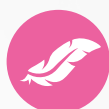

EXTRAORDINARILY LIGHT

The device weighs only 99 g, which makes it easy to use for treatments.



**FROM ODONTOLOGISTS
TO ODONTOLOGISTS**

TORQUE CONTROL

In different working modes the **torque value can be adjusted differently.**

Z-Evo mini manages file stress in a smart way.

When the torque limit is reached, instead of going into reverse mode, the motor oscillates until the stress is reduced. Once the torque level returns to the predetermined setting, the file goes back to continuous rotation.

In M0 reciprocating mode, **5 pre-set reciprocating angle values** can be accessed.



TORQUE CALIBRATION
For an accurate fit.



OLED SCREEN
For a higher-quality image.



AMBIDEXTROUS MONITOR
The screen can be rotated to be used with the left or right hand.



THE LIGHTEST ENDODONTIC MOTOR ON THE MARKET

It weighs less than 100 g, reducing hand fatigue of the endodontist during treatment and providing an extraordinary user experience.

BETTER VISIBILITY AND REAR ACCESS

Its rotary contra-angle is small in size.



Z-Evo mini contra-angle with reduced-size head and neck.



Conventional contra-angle.

ZEvo mini

WELCOME TO ENDODONTICS OF THE FUTURE

We are universalizing endodontics,
promoting accessibility to a safe
and quality working environment.

Technical Specifications

Speed Range	120 - 650 rpm
Torque Range	0.5 - 4.0 N•cm
Number of memories	10 (M0 - M9)
Battery Power	800 mAh



Z-Evo mini Reference Codes

Z-Evo mini Endodontic Motor **6080240**

Z-Evo mini Contra-angle **6036101**



FROM ODONTOLOGISTS
TO ODONTOLOGISTS

hola@zarc4endo.com
www.zarc4endo.com

