




Pluraject™

- en Ampoule Syringe
- de Ampullenspritze
- fr Seringue à cartouche
- it Siringa per cartucce
- es Jeringa para ampollas
- pt Seringa para cartuchos
- nl Ampulinjectiespuit

rev. 00 - 2024/09



Pluraject™ 2

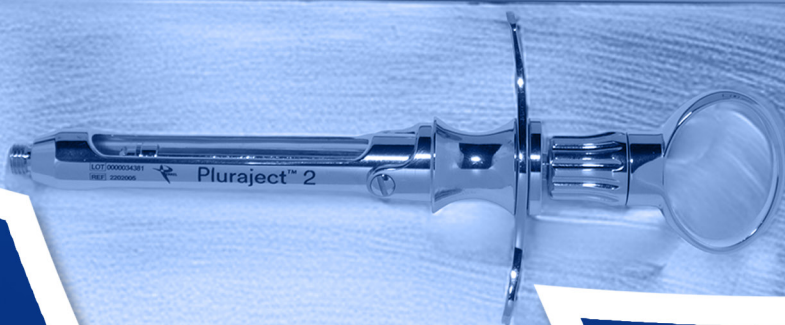

Italian quality and
innovation for oral health

PLURAJECT 2™

Aspiration syringe

Properties

- › Stainless steel syringe rack
- › Syringe rack with flip mechanism; the ampoule holder is flipped aside, and only the ampoule itself is exchanged
- › M6 thread (European market) and 40 TPI (British market)
- › Aspiration
 - Manual aspiration by means of three extending aspiration hooks
 - Self aspiration
- › Reprocessing (cleaning, disinfection, and sterilization) has been validated for 500 times.

Pluraject 2 Syringe – Cartridge Engagement System:

- The Pluraject 2 syringe is designed for use with dental anesthetic cartridges featuring a rear elastomer plunger with a central cavity (designed for harpoon engagement).
- Compatible with dental cartridges with perforable/engageable plungers (e.g., Ubistesin).



Advantages

- › Compact and lightweight for better handling
- › Cartridges are easy to insert and remove thanks to the flip mechanism
- › Appropriate for both active and passive aspiration
- › Inch and metric thread available



Italian quality and innovation for oral health

Founded in 1948, Pierrel is an Italian pharmaceutical company based in Capua, specializing in the development, production, registration, and licensing of dental anesthetics and medical devices for oral health.

www.pierrelgroup.com

BINO PHOTO