

## Simple, reliable extraction

Typenex® Surgical offers stone extractors and antimigration devices to aid in the removal of kidney stones for use in both ureteroscopy and percutaneous procedures.

**TYPENEX**   
SURGICAL

# Stone Extractors



## RocaStone Nitinol Tipless Basket


Available in 1.3 Fr, 1.9 Fr, and 2.2 Fr, the RocaStone Tipless baskets are ideal for the removal of upper urinary tract and renal stones. The small size of the basket improves irrigation and visibility.

The tipless configuration allows for close stone proximity to successfully retrieve calyceal calculi. The distal surface is flat, eliminating tissue to tip interface and resulting in less traumatic manipulation. The flexible shaft allows for enhanced scope deflection.



## RocaStone Flat & RocaStone Helical

Available in 2.5 Fr and 4.0 Fr, the RocaStone 4 Wire baskets are designed for rigid and semi-rigid ureteroscopy procedures. Both the flat wire and the helical design feature ergonomic, removable handles. Choose the wire design that best meets your stone manipulation needs.

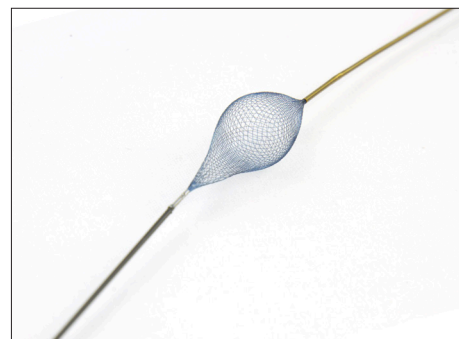

## XenX™

Developed to trap fragments larger than 1 mm, the XenX™ immobilizes stones for further disintegration and/or removal by basket or irrigation. Available in 0.038 in, its use aims to shorten procedure time, reduce invasiveness, and address the clear demand for a cost effective stone retention device.

## Ordering Information

Product Code	Description
ROSH2500ST	RocaStone “Helical” Nitinol Stone Extractor 2.5 Fr x 120 cm, 16 mm basket
ROSH4000ST	RocaStone “Helical” Nitinol Stone Extractor 4.0 Fr x 120 cm, 16 mm basket
ROSS2500ST	RocaStone “Flat” Nitinol Stone Extractor 2.5 Fr x 120 cm, 16 mm basket
ROSS4000ST	RocaStone “Flat” Nitinol Stone Extractor 4.0 Fr x 120 cm, 16 mm basket
ROSZ1300ST	RocaStone Nitinol Stone Extractor Tipless 1.3 Fr x 120 cm, 12 mm basket
ROSZ1900ST	RocaStone Nitinol Stone Extractor Tipless 1.9 Fr x 120 cm, 16 mm basket
ROSZ1920ST	RocaStone Nitinol Stone Extractor Tipless 1.9 Fr x 120 cm, 12 mm basket
ROSZ2200ST	RocaStone Nitinol Stone Extractor Tipless 2.2 Fr x 120 cm, 16 mm basket
ROXE03812ST	XenX™ – Stone Management Device (12 cm tip length, mid-distal stone), 150 cm length
ROXE03807ST	XenX™ – Stone Management Device (7 cm tip length, proximal stone), 150 cm length

For more information, visit [www.typenexsurgical.com](http://www.typenexsurgical.com).

Products distributed by Typenex® Surgical. RocaStone products manufactured by Promepla. XenX™ manufactured by Xenolith.