


Access redefined

Designed for gaining access to a tortuous ureter, the RioTracer® Wire Guide has a stiff nitinol core and a hydrophilic, barium-filled, radiopaque, floppy tip to enhance fluoroscopic visualization.

RioTracer®

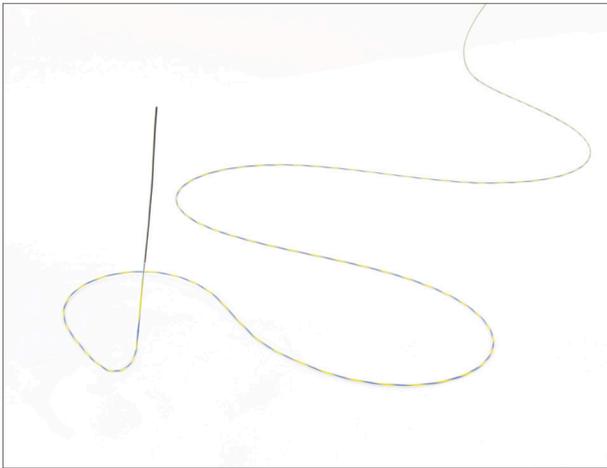

Hydrophilic Floppy Tip – Facilitates advancement beyond obstructions and negotiates challenging anatomy.

Nitinol Core – Provides resistance to kinking for enhanced control.

PTFE Coating – Allows for improved handling and facilitates enhanced catheter placement.

Striped Pattern – Designed to provide clear endoscopic visualization of wire movement.

Barium-filled Radiopaque Tip – Enhances visualization under fluoroscopy.

Ordering Information

Product Code	Description
ROWR3500ST	RioTracer® Hydrophilic Dual Flexible, 0.035 in x 150 cm, 3 cm tip
ROWR3501ST	RioTracer® Hydrophilic Angled/Straight Dual Flexible, 0.035 in x 150 cm, 3 cm tip
ROWR3502ST	RioTracer® Hydrophilic and Floppy Tip, 0.035 in x 150 cm, 3 cm tip

For more information, visit www.typenexsurgical.com.
Products distributed by Typenex® Surgical. Manufactured by Promepla.