


2015-04-30

OVUCHECK® Reader

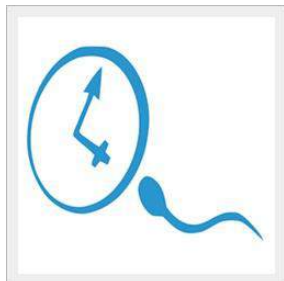

PRODUCT CODE: TRD-751

DESCRIPTION

OVUCHECK® Reader is a Convertible Well-strip/Cuvette Reader and is a portable miniphotometer designed to do ELISA testings in either 1x8/1x12 microwell strips or round/square cuvettes by interference filter for wavelength of 400-700nm.

OVUCHECK® Reader is a user-friendly ELISA reader. All operations are carried out with only two touch buttons by following instruction on the LCD display, and results can be obtained within two seconds.

OVUCHECK® Reader undertakes testings with an extreme economic Reagent consumption, in microwell at 0.3ml and in cuvette at 0.5ml.

SUMMARY AND EXPLANATION

OVUCHECK READER® can be use with Biovet Ovucheck® kits to measure progesterone levels in milk or blood.

CHARACTERISTICS

- portable, compact, accurate, easy to use and affordable
- Holders for 12mm diameter tubes or 9 x 13mm well strips.
- Houses multiple filters, including the 405 nm and 550 nm filters
- Changing filters and holders is fast and very easy.
- Instructions that are as simple as using a video device.
- Measurement accuracy comparable to more expensive devices.



Biovet inc.
4375, Beaudry Ave
St-Hyacinthe QC J2S 8W2
Canada
T 1-888-8BIOVET
F 450 771-2996
order@biovet-inc.com

www.biovet.ca