



ALLIANCE  
BIO EXPERTISE

# MEDIAWEL 30<sup>®</sup>

Media Preparator





# MEDI AWEL30®

Increase your productivity  
quality  
traceability

## PERFORMANCES

- 2.5 to 30 Liters
- 30 – 125 °C  
+/- 0.5 °C
- 50 – 90 mins

CONNECTED  
MICROBIOLOGY

# MEDI AWEL®

Safety lid  
(locked above 80°C)




Uphill & downhill  
Filtration  
Contamination prevention

Double pressure monitoring



Large stirrer  
(wide viscosity range)

Medical wheels

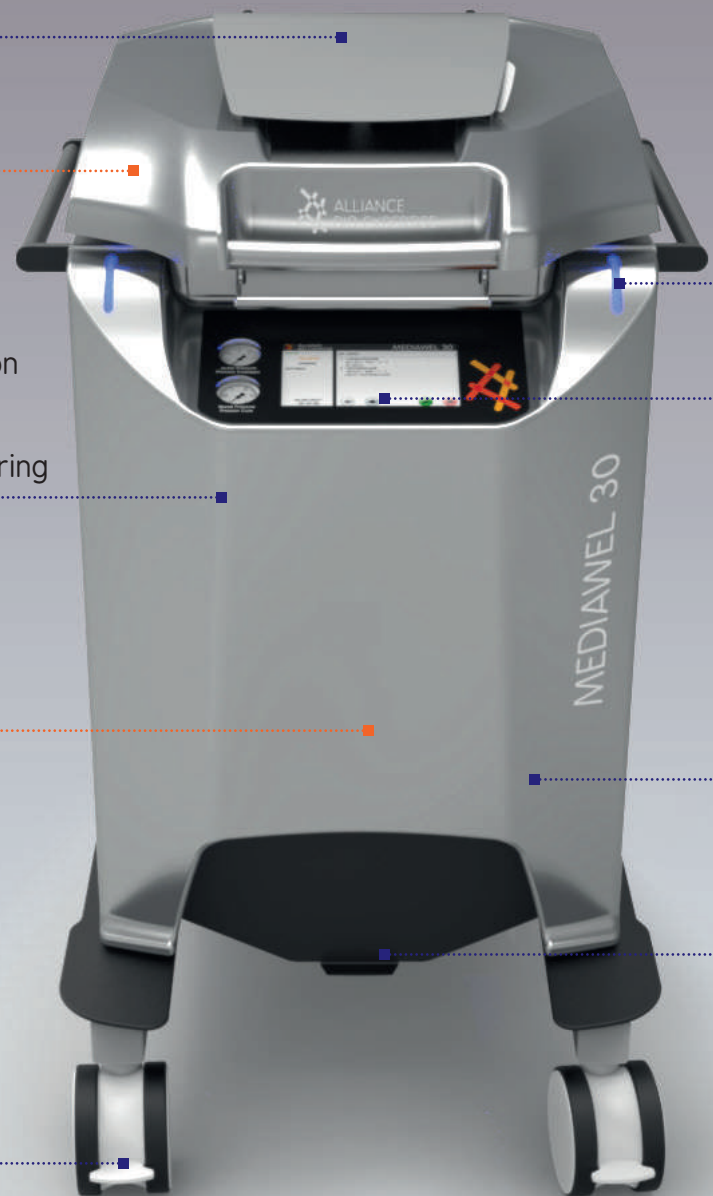

WIFI live communication

Multicolour visual alarms

7" touch screen LCD  
Intuitive display



Anti-shock body  
(ABS PC PLASTICS)

Trapdoor  
(easy cleaning of the bowl)



# MEDI AWEL30®

Increase your productivity  
quality  
traceability

ALLIANCE BIO EXPERTISE  
WIDE RANGE OF MEDIA



Ask for the list!  
[www.alliance-bio-expertise.com](http://www.alliance-bio-expertise.com)

CONNECTED  
MICROBIOLOGY

## SIMPLIFY



ALL IN ONE  
MEDIA PREPARATION

## QUALITY MEDIA

- Controlled temperature = perfect sterility – Perfect Fertility  
→ T° precision: +/- 0.5°C
- Continuous homogenization
- Large smooth stirrer (1807 cm<sup>2</sup>) for higher homogeneity. Slow stirring (9 adjustable speeds 20 – 110 rpm) to avoid bubbling
- Less handling = less contamination

## TIME SAVING

Before



With MEDI AWEL




## TRACEABILITY



Temperature monitoring  
Up to 3 points

eTRACE®

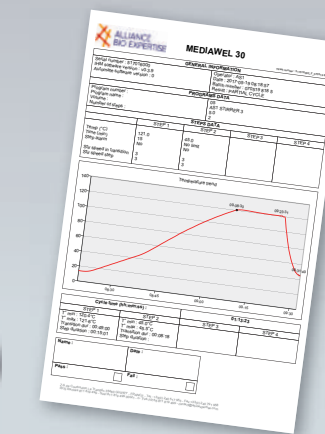







## ALLIANCE BIO EXPERTISE



Empowered from our people's skills and experiences, we, ALLIANCE BIO EXPERTISE, together with our customers and partners decided to manufacture Microbiology automated equipment around essential values:  
Customer Focus, Service, Loyalty, Innovations

**Ask your local partner  
for a personal demonstration!**

ZI de courbouton , Le Tremplin – 35480 GUIPRY - FRANCE  
E-mail : [contact@abioexpertise.com](mailto:contact@abioexpertise.com)  
Phone : +33(0)2.40.51.79.53  
[alliance-bio-expertise.com](http://alliance-bio-expertise.com)