

JumpStart® Bunny®

SMO trimline with proximal PF resist strap.

CASCADE®
daFO®

The JumpStart Bunny was developed to help young patients with fully-correctable pronation, supination, or swing phase inconsistency to realign their base of support for improved stability in standing and walking. Full wrap-around flexible support encompasses the foot and ankle for pronation and supination control. A well-defined plantar surface and wrap-around compression provide proprioceptive input. Proximal strap(s) provides an aid to frontal plane control and influence mild resistance to excess sagittal plane motion. Solid outer foundation adds extra arch support and heel stability.



Key Features

High proximal "ears"

provide additional medial/lateral control.

Proximal posterior elastic strap

adds to frontal plane control and offers mild resistance to excess plantarflexion. Option to replace posterior strap with anterior/posterior wrap-around strap.

Solid outer foundation

provides additional support and heel stabilization.

Full wrap inner liner

made from soft, thin polyethylene plastic allows full alignment control of the heel, mid-foot, and forefoot.

Choice of two outer foundation materials

Polyethylene for moderate flexibility or co-poly for firm support.

Prefabricated

sized by measurement using the Fast Fit sizing jig.

Widths: Wide and Narrow

Sizes: 4.00 – 9 inches (0.25 inch increments)

Sizes in Össur Nordic stock: 4.50-7.75 in. wide

Designed for patients who:

- Exhibit a fully correctable foot position and fit within the product size range.
- Present with moderate to strong pronation or supination and associated gait instability.
- May benefit from the sensory feedback provided by a fully contoured foot bed.
- Demonstrate mild swing phase instability, or excess plantarflexion presented as mild tow walking or knee hyperextension.
- Require some resistance to plantarflexion.



Wrap-around strap

Co-poly



No casting needed!

  
WWW.OSSUR.SE

DISTRIBUTED BY

Össur Nordic
Box 67
751 03 Uppsala
Sweden

TEL +46 (0)18 18 22 00
FAX +46 (0)18 18 22 18
info@ossur.com

 **ÖSSUR**®
LIFE WITHOUT LIMITATIONS




No Casting

High SMO with proximal strap, PF resist, optional DF resist

Today's Date: _____

Patient	Last name:			
	First name:	<input type="checkbox"/> Male <input type="checkbox"/> Female		
	Birth date:			
	Parent or Guardian:			
Practitioner	Name:	Title:		
	Facility:			
	Street address:			
	City:	State:	Zip:	
	Phone:			
	Email:			
Payment Options	<input type="checkbox"/> Facility Billing (Practitioner)	-OR-		
	Account Name or #:			
	P.O. N°:	<input type="checkbox"/> CC on file		
	<input type="checkbox"/> Insurance Billing (Parent / Guardian / Practitioner)	-OR-		
	—UCAN N°:			
	<input type="checkbox"/> Direct Purchase (Parent / Guardian)			
	<input type="checkbox"/> Check attached			
	Credit Card: <input type="checkbox"/> Visa <input type="checkbox"/> MasterCard <input type="checkbox"/> AMEX <input type="checkbox"/> Discover			
	Cardholder's Phone:			
	Credit Card No:			
Billing Information	Billing Name:			
	Facility:			
	Street address:			
	City:	State:	Zip:	
	Phone:			
Shipping	<input type="checkbox"/> Same as billing information.	-OR-		
	Shipping contact name:			
	Street address:			
	City:	State:	Zip:	
Phone:				

Size | Outer Shell | Straps | Options

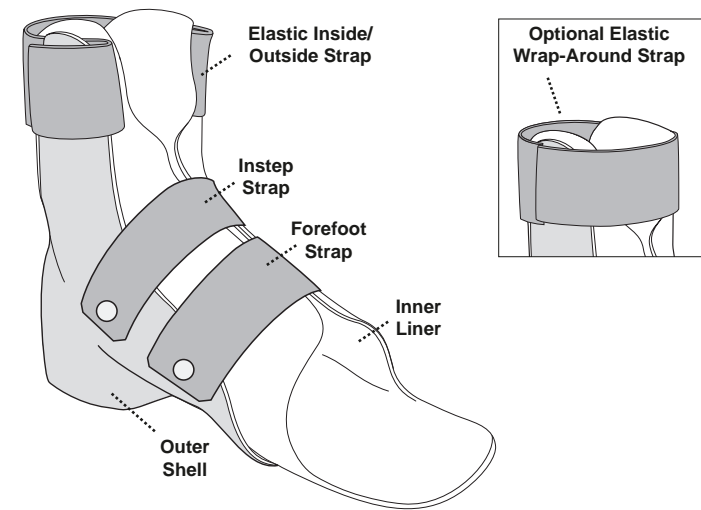

1 Sizing Pair Left Right

Length: _____ 4.00 – 9.00 in. (0.25 in. increments)

Width: Wide Narrow

2 Outer Shell **Moderate Flexibility – Polyethylene**
Recommended for sizes 4.00 – 5.75 (available for all sizes)

Shell color: Blue Pink

-OR- **Firm – Co-poly** (shell color: **White only**)
Recommended for sizes 6.00 – 9.00 (available for all sizes)

3 Straps

Color: Blue Pink

Posterior: choose one Elastic Inside / Outside (resists PF)
 Elastic Wrap-Around (resists PF and DF)

Instep: choose one Riveted layover
 Layover (no rivets)
 Riveted D-ring

Forefoot: choose one Riveted Layover
 Layover (no rivets)

4 Options Instep pad
 Toe rise pad
 Toe rise pad with abduction strap

Comments

© 2015 Cascade Dafo, Inc. All rights reserved.

Thank you!