

Fast Fit[®] Cricket[®]



CASCADE[®]
daFO[®]

Sub-malleolar, UCBL trimline

The Cricket was developed especially for children with mild to moderate pronation and supination. Stability and alignment improvement come from a deep heel cup, well-defined plantar surface, and medial / lateral forefoot containment in combination with a stable outer foundation. The UCBL trimline allows the orthosis to become enclosed entirely within a shoe.



Key Features

Solid outer foundation

provides additional arch support and heel stabilization.

Partial wrap inner liner

made from soft, thin, polyethylene plastic cups heel, supports arches, and contains mid-foot and forefoot. Trimline is below top edge of shoe.

Choice of two outer foundation materials

Polyethylene for moderate flexibility or co-poly for firm support.

Prefabricated

sized by measurement using the Fast Fit sizing jig.

Widths: Wide and Narrow

Sizes: 4.00 – 9 inches (0.25 inch increments)

Sizes in Össur Nordic stock: 4.50-7.75 in. wide

Designed for patients who:

- Exhibit a fully correctable foot position and fit within the product size range.
- Present with mild to moderate pronation or supination and associated gait instability.
- May benefit from the sensory feedback provided by a fully contoured foot bed.
- Require access to free dorsiflexion and planterflexion.



No casting needed!



WWW.OSSUR.SE

DISTRIBUTED BY

Össur Nordic
Box 67
751 03 Uppsala
Sweden

TEL +46 (0)18 18 22 00
FAX +46 (0)18 18 22 18
info@ossur.com



ÖSSUR[®]
LIFE WITHOUT LIMITATIONS




No Casting

Today's Date: _____

Patient	Last name:		
	First name:	<input type="checkbox"/> Male	<input type="checkbox"/> Female
	Birth date:		
	Parent or Guardian:		
Practitioner	Name:	Title:	
	Facility:		
	Street address:		
	City:	State:	Zip:
	Phone:		
	Email:		
Payment Options	<input type="checkbox"/> Facility Billing (Practitioner)	-OR-	
	Account Name or #:		
	P.O. N°:	<input type="checkbox"/> CC on file	
	<input type="checkbox"/> Insurance Billing (Parent / Guardian / Practitioner)	-OR-	
	—UCAN N°:		
	<input type="checkbox"/> Direct Purchase (Parent / Guardian)		
	<input type="checkbox"/> Check attached		
	Credit Card: <input type="checkbox"/> Visa <input type="checkbox"/> MasterCard <input type="checkbox"/> AMEX <input type="checkbox"/> Discover		
	Cardholder's Phone:		
	Credit Card No:		
Billing Information	Exact name on card:		
	Exp. Date:	V-code:	
	Billing Name:		
	Facility:		
	Street address:		
	City:	State:	Zip:
Shipping	Phone:		
	Email:		
	<input type="checkbox"/> Same as billing information.	-OR-	
	Shipping contact name:		
Street address:			
City:	State:	Zip:	
Phone:			

Size | Outer Shell | Options

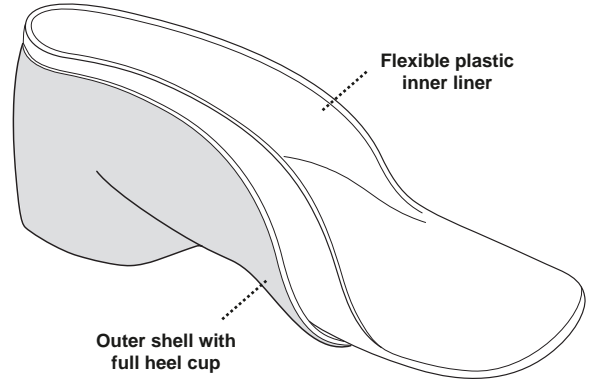

1 Sizing Pair Left Right

Length: _____ 4.00 – 9.00 in.
(0.25 in. increments)

Width: Wide Narrow

2 Outer Shell **Moderate Flexibility – Polyethylene**
Recommended for sizes 4.00 – 8.00 (available for all sizes)

Shell color: Blue Pink

-OR-

Firm – Co-poly (shell color: White only)
Recommended for sizes 8.25 – 9.00 (available for all sizes)

3 Options Toe rise pad

Toe rise pad with abduction strap

Comments

Order Cricket Rev.10 (Oct 2013)

Thank you!